

URIS OMNI

URIS IMPLANTS

IMPLANTS

URIS IMPLANTS

Powered by **TruAbutment**



WHAT IS
URIS



ULTIMATE
RESTORATIVE
IMPLANT
SOLUTION



ULTIMATE
RESTORATIVE
IMPLANT
SOLUTION




URIS Prosthetics & Implant Kit

1.2.3. DONE!

1

Clinician sends captured data

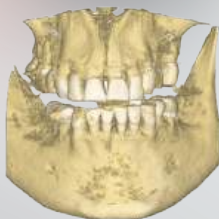

☒ CT Scan



☒ Intraoral Scan (IOS)

& orders **UPIK**
on shop.truabutment.com

2

URIS manufactures:



Surgical Guide



Prosthetics

3



Day of Surgery

Clinician is ready for guided surgery
& delivery of temporary prosthetics





PICK!

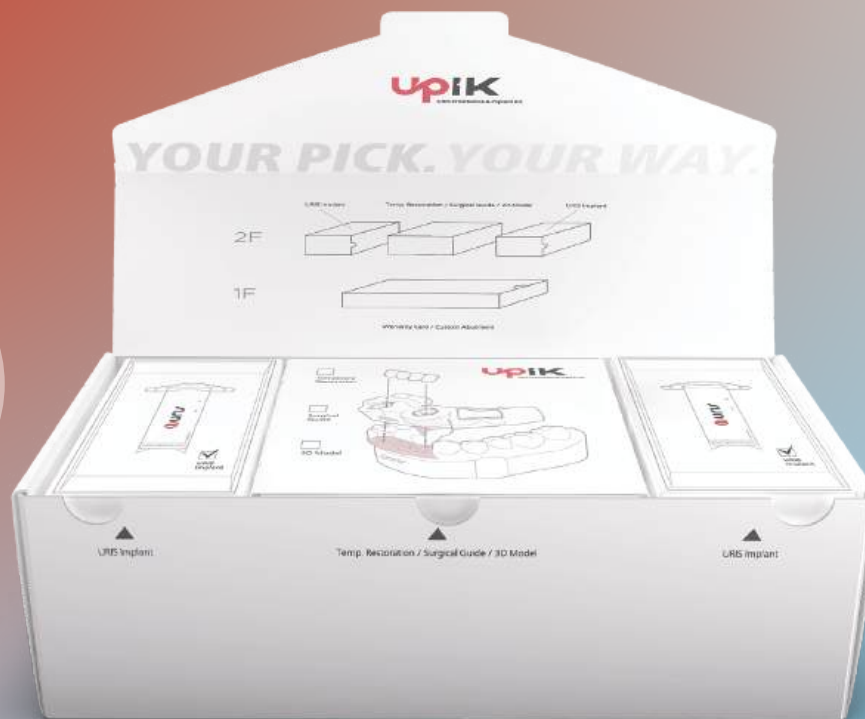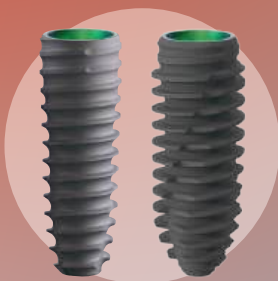

Start by selecting your choice of **implant & prosthesis**. We will put together everything you need from **start ► finish** in a box, ready for you to deliver chairside.

NEW

Open System 

Compatible with any dental CAD software for data importation and design.

Minimized Scanning Area 

Focus on scanning **just the blue area** to reduce scanning error. Easy access and shorter scanning time.



ioConnect™

 Effortless MUA Detection with your IOS

Empower your **existing IOS (Intra-Oral Scanner)** to achieve the same photogrammetry workflow.

 All About Options

Multiple sizes (S, M, L) available for you to choose depending on the implant position and the anatomy of each patient.

Effortless MUA Detection with Your IOS

ioConnect™ KIT +  **SPECIAL OFFER**
GIFT 1 free pack of Direct Screw (20 screws)

URIS (M1.6)

MSRP \$3,000




HS caps x 12 Short/Medium/Long MUA Detector x 4 x 4 x 4 1.27 hex hand driver x 1

 **Direct Screw** x 20



Components


Item Code	Description
HS-UR	HS Cap compatible with URIS
HSS-UR	HS Cap Small compatible with URIS
IO-URL	ioConnect compatible with URIS Long
IO-URM	ioConnect compatible with URIS Medium
IO-URS	ioConnect compatible with URIS Short
IOK-UR	ioConnect Kit compatible with URIS

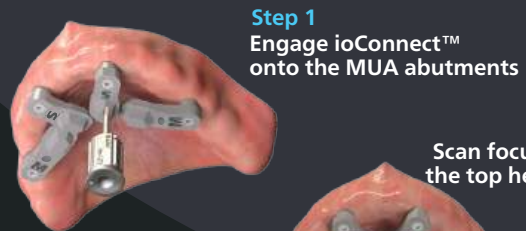

Step 1
Engage ioConnect™
onto the MUA abutments



Step 2
Scan focusing on
the top hexagons

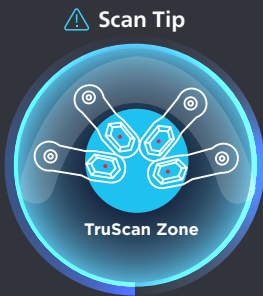

⚠ Scan Tip

TruScan Zone



Step 3
Choose your preferred HS Cap
(2 types available)



Small cap is included in the full kit



Sold separately



Step 4
Scan HS Caps to capture as
soft tissue reference points

⚠ Scan Tip >>> Focusing specifically on the hexagonal portion
and **minimizing the scanning area** will improve accuracy.

Do



Do Not



TruSuite

trusuite.truabutment.com




ioConnect™



Scan to **learn** more

After scanning the ioConnect™ scan bodies, export it and
import it to the **ioConnect™ app**. Then you can import the
HS cap scan (gingiva data) and the exported ioConnect™ (MUA data)
to any dental CAD software to design your temporary/final arch.

THIS IS IT. WE ARE TRU.

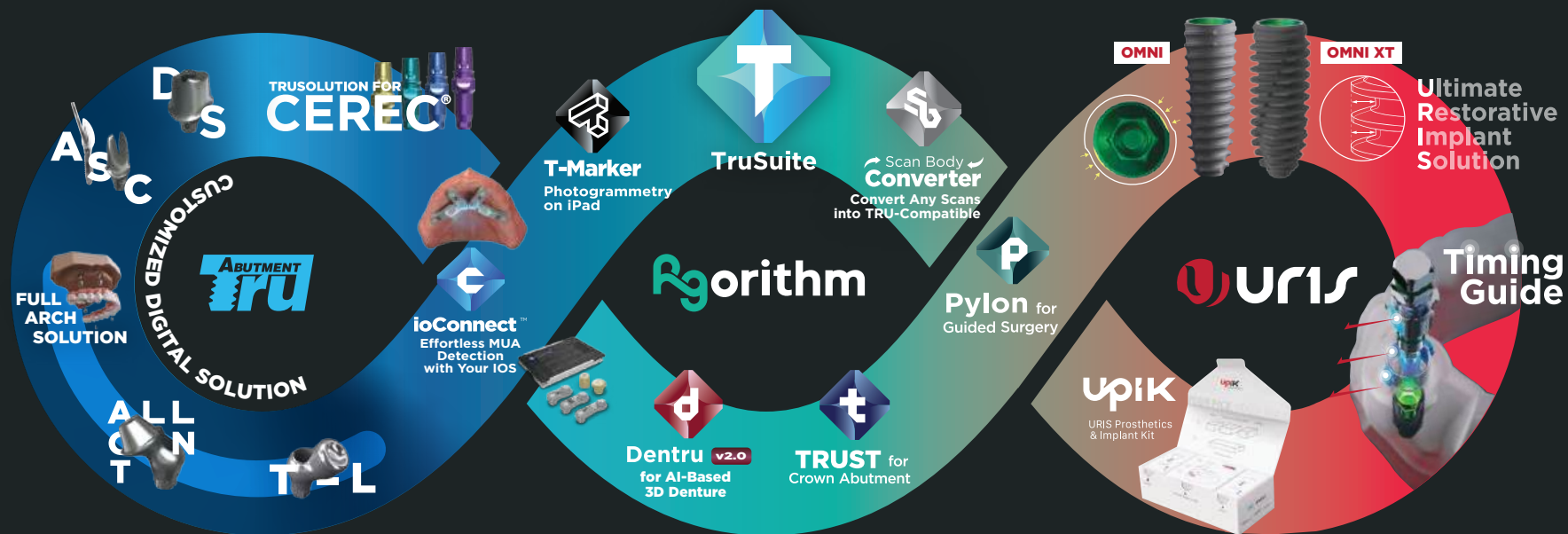

PYLON
for Guided Surgery

Dentru
for 3D Denture

MBuilder
for 3D Model

NGuard
for Nightguard

Scan Body
Converter
Convert Any Scans
into TRU-Compatible

NEW

ioConnect™
Effortless MUA Detection
with Your IOS

HOT

T-Marker
Photogrammetry
on iPad

COMING SOON

TRUST
for Crown & Abutment



TruSuite
Level Up Now.
Be a part of the new big thing.

CONTENTS



FIXTURE

- 1. OMNI PG.22
- 2. OMNI Tapered PG.24
- 3. Cover Screw PG.24



CAD/CAM SOLUTION

- 1. URIS Scan Body PG.28
- 2. URIS Base PG.28
- 3. Digital Lab Analog PG.28
- 4. Multi-Unit Digital Components PG.29
- 5. Surgical Guide Sleeve PG.29




URIS SOLUTION FOR CEREC®

- 1. URIS Scan Post for CEREC® PG.32
- 2. URIS Base for CEREC® PG.32
- 3. Alpha and Beta Kit PG.33



PROSTHETICS

- 1. Healing Abutment **PG.36**
- 2. Temporary Abutment **PG.39**
- 3. Direct Abutment **PG.41**
- 4. Basic Abutment **PG.45**
- 5. Milling Abutment **PG.52**
- 6. Angled Abutment **PG.55**
- 7. Multi-Unit Abutment **PG.58**
- 8. T-L Abutment **PG.61**
- 9. Ball Abutment **PG.63**

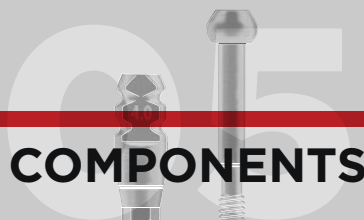

COMPONENTS

- 1. Impression Components **PG.68**
 - 2. Protective Cap **PG.72**
 - 3. Impression Cap **PG.72**
- 4. Abutment Components **PG.72**
- 5. Multi-Unit Components **PG.73**
- 6. T-L Components **PG.75**
- 7. Ball Components **PG.75**
- 8. Fixture Analog **PG.75**
- 9. Direct Abutment Analog **PG.76**
- 10. Ball Abutment Analog **PG.77**
- 11. Abutment Screws **PG.77**



INSTRUMENTS


- 1. Surgical Kit **PG.80**
- 2. Wide Surgical Kit **PG.92**
- 3. Tapered Surgical Kit **PG.98**
- 4. Surgical Plus Kit **PG.102**
 - 5. Pylon Kit **PG.106**
 - 6. Pylon Plus Kit **PG.118**
- 7. Tapered Pylon Kit **PG.122**
- 8. Narrow Pylon Kit **PG.126**
- 9. Pylon Anchor Kit **PG.134**
- 1. Pylon Crestal Sinus Kit **PG.138**
 - 2. Prosthetic Kit **PG.148**
 - 3. ASC Kit **PG.152**

ULTIMATE RESTORATIVE IMPLANT SOLUTION


FLOWCHARTS



Fixture Level Impression

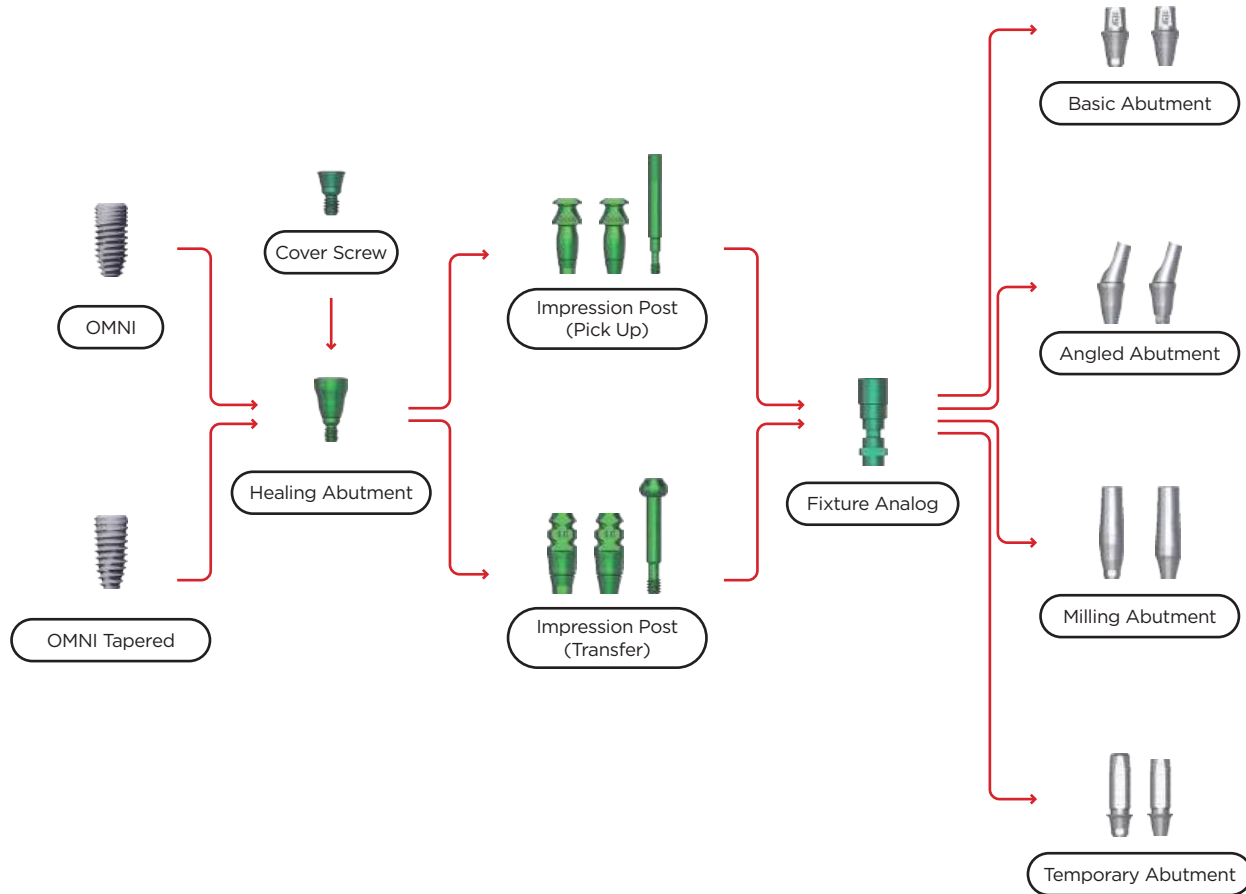


Digital Impression & Components

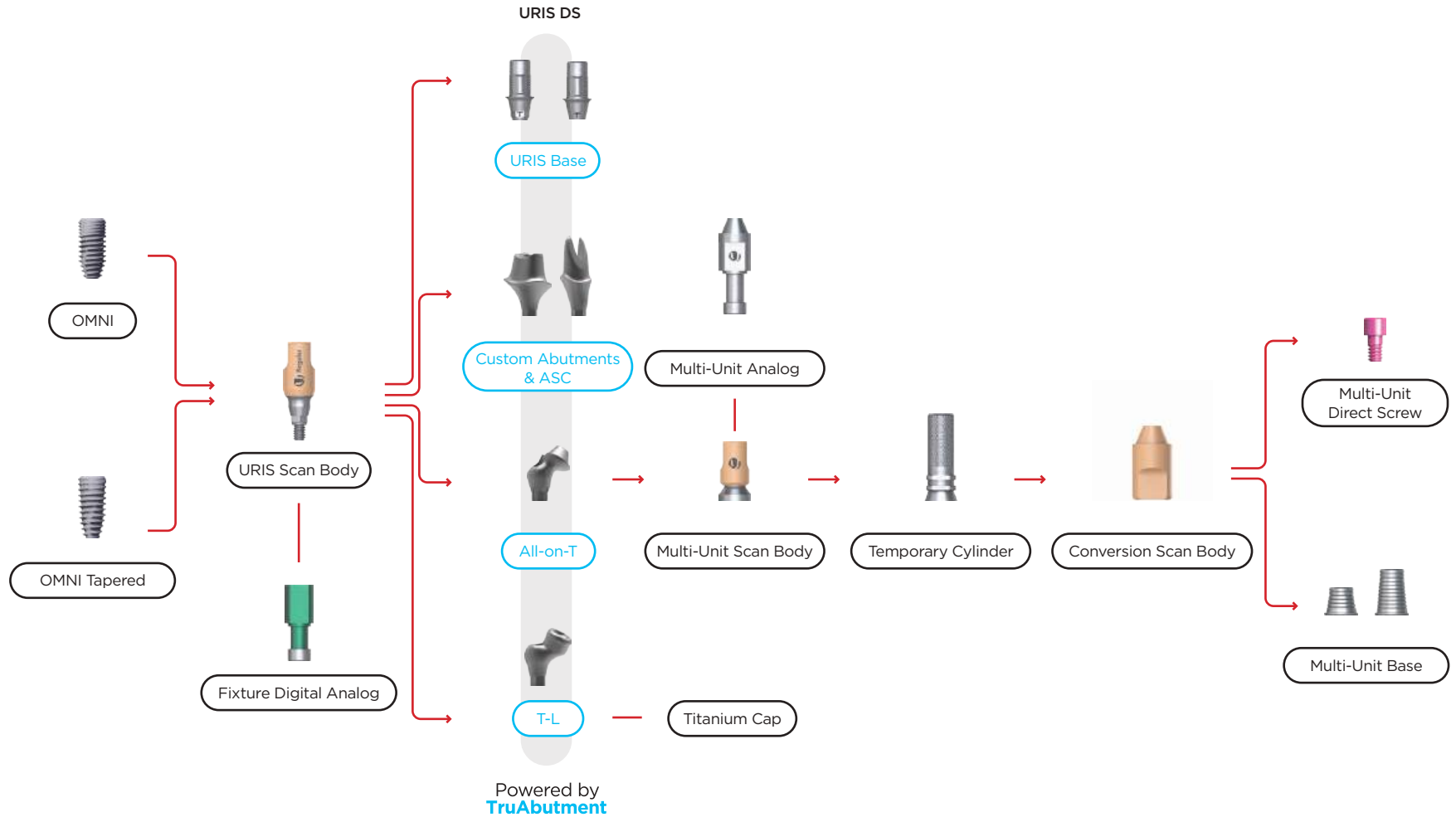


Abutment Level Impression

Overdenture - T-L Abutment

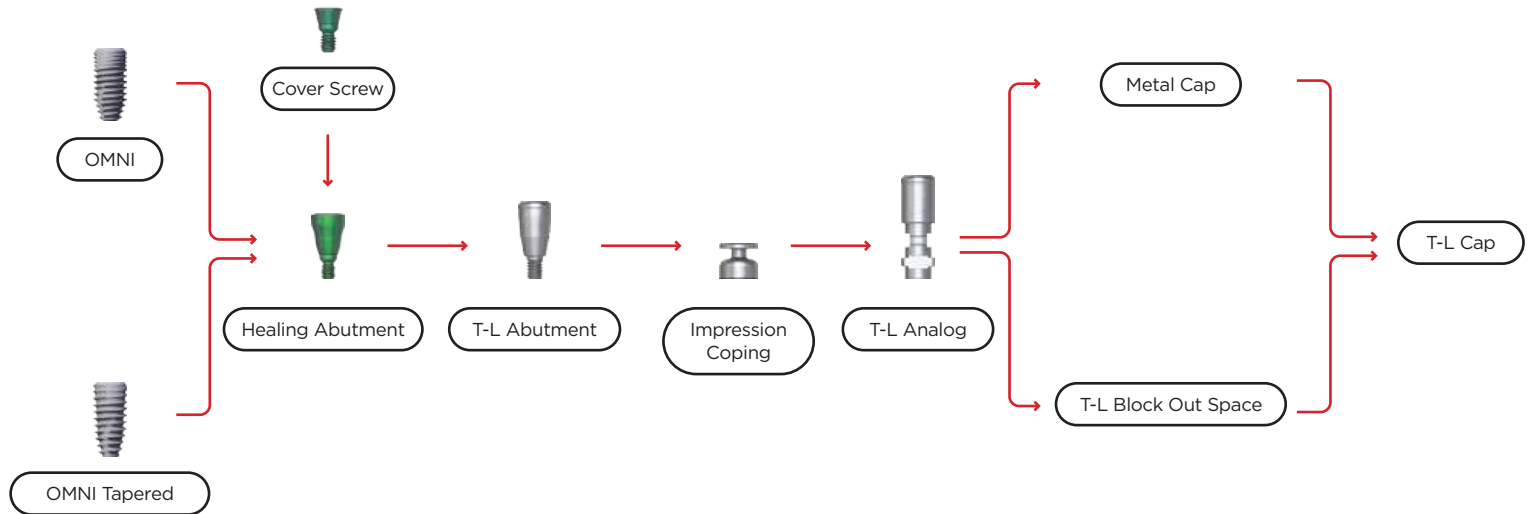

Abutment Level Impression

Multi-Unit Abutment & Components

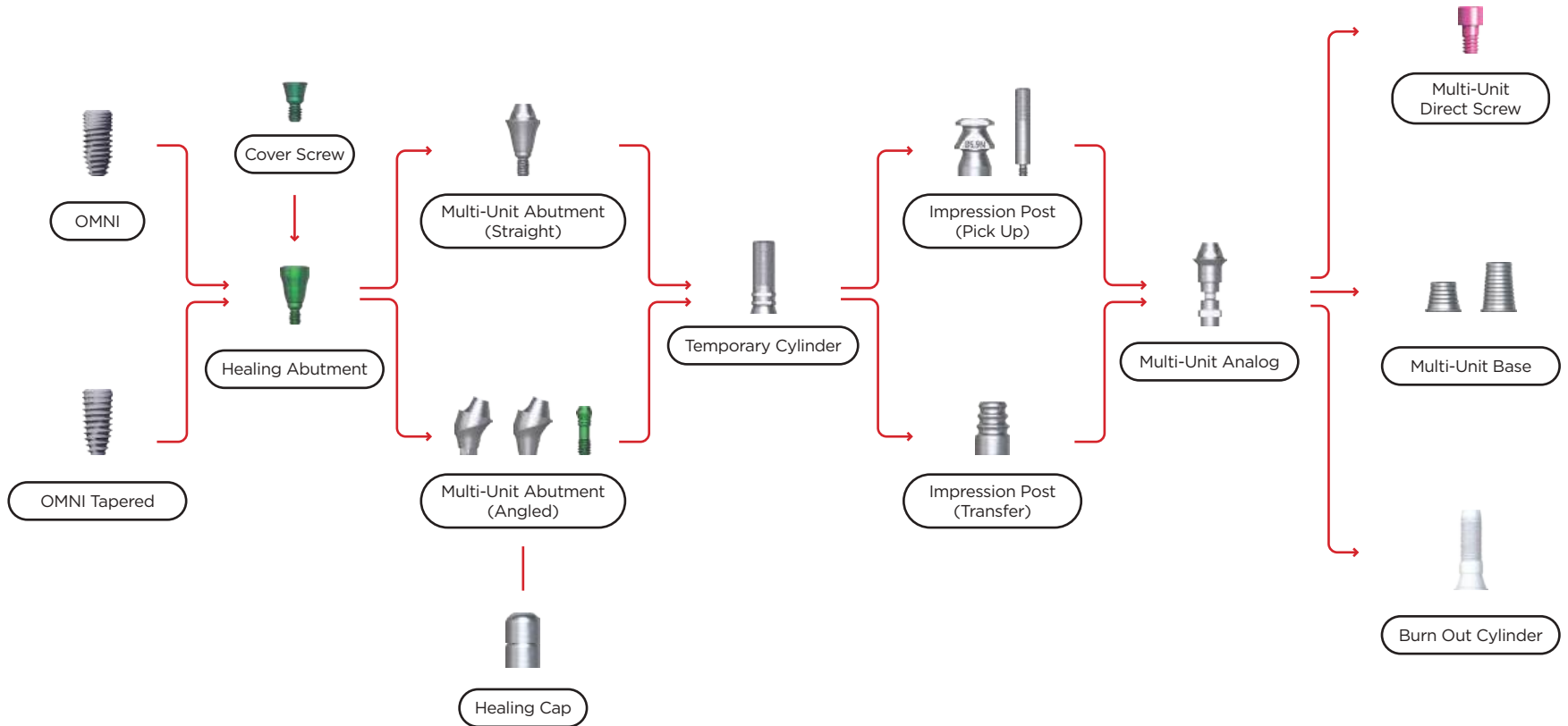

FIXTURE

CAD/CAM SOLUTION


URIS SOLUTION FOR CEREC®

PROSTHETICS

COMPONENTS

INSTRUMENTS

1. OMNI
2. OMNI Tapered
3. Cover Screw

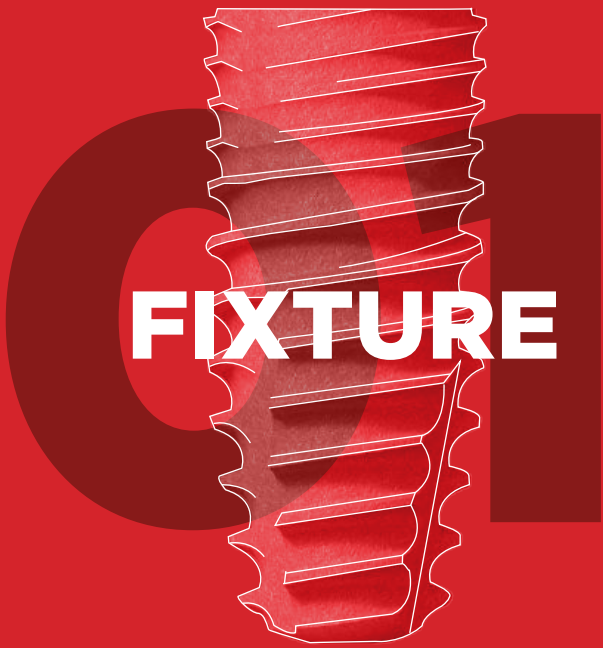

OMNI

***1 Cover Screw Included**

Platform	Diameter								
1.95 Narrow	D Ø 3.0	10.0	11.5	13.0	14.5				
	Item No.	OMF30 100	OMF30 115	OMF30 130	OMF30 145				
D Ø 3.5	8.5	10.0	11.5	13.0	14.5	16	18		
	Item No.	OMF35 085	OMF35 100	OMF35 115	OMF35 130	OMF35 145	OMF35 160	OMF35 180	
2.5 Regular	7.0	8.5	10.0	11.5	13.0	14.5	16	18	
	Item No.	OMF40 070	OMF40 085	OMF40 100	OMF40 115	OMF40 130	OMF40 145	OMF40 160	OMF40 180
D Ø 4.5	7.0	8.5	10.0	11.5	13.0	14.5	16	18	
	Item No.	OMF45 070	OMF45 085	OMF45 100	OMF45 115	OMF45 130	OMF45 145	OMF45 160	OMF45 180

OMNI

*1 Cover Screw Included

Scale 150%

Unit : mm



Regular

Platform	Diameter						
 D Ø 5.0		 7.0	 8.5	 10.0	 11.5	 13.0	 14.5
	Item No.	OMF50070	OMF50085	OMF50100	OMF50115	OMF50130	OMF50145
 D Ø 5.5		 7.0	 8.5	 10.0	 11.5		
	Item No.	OMF55070	OMF55085	OMF55100	OMF55115		
 D Ø 6.0		 7.0	 8.5	 10.0			
	Item No.	OMF60070	OMF60085	OMF60100			
 D Ø 6.5		 7.0	 8.5	 10.0			
	Item No.	OMF65070	OMF65085	OMF65100			

01 Fixture

02 CAD/CAM Solution

03 URIS Solution for CEREC®

04 Prosthetics

05 Components










06 Instruments

23

URIS IMPLANTS

OMNI Tapered

*1 Cover Screw Included

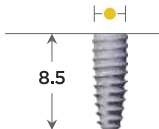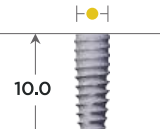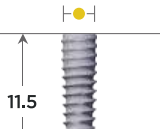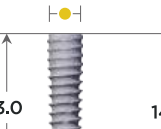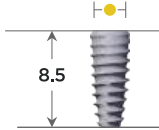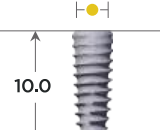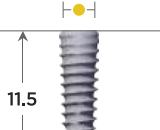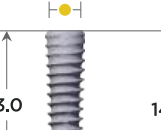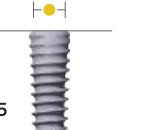Platform	Diameter					
Scale 150% Unit : mm  Narrow	 D Ø 3.0	Item No.	 8.5 NTF30 085	 10.0 NTF30 100	 11.5 NTF30 115	 13.0 NTF30 130
	 D Ø 3.5	Item No.	 8.5 NTF35 085	 10.0 NTF35 100	 11.5 NTF35 115	 13.0 NTF35 130
			 14.5 NTF35 145			

Cover Screw

Hex driver: 1.27 | Torque: 5-10 N·cm

Platform	
N Narrow	 UCS16
R Regular	 UCS20



















OMNI Tapered

*1 Cover Screw Included

Scale 150%
Unit : mm



Regular

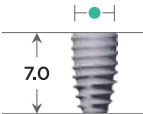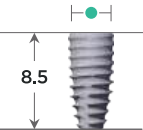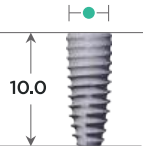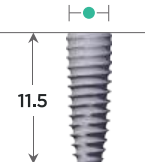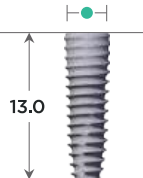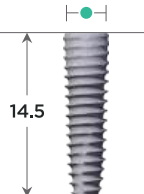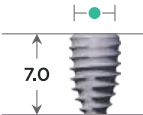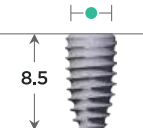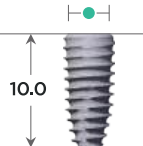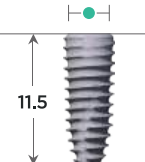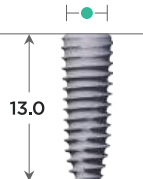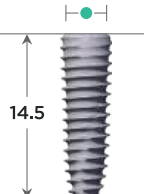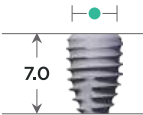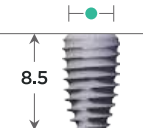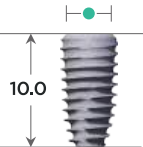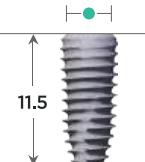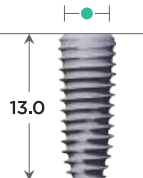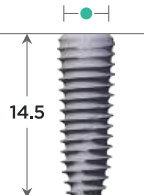Platform	Diameter						
 Regular	 D Ø 4.0	 7.0	 8.5	 10.0	 11.5	 13.0	 14.5
	Item No.	NTF40 070	NTF40 085	NTF40 100	NTF40 115	NTF40 130	NTF40 145
 D Ø 4.5	 7.0	 8.5	 10.0	 11.5	 13.0	 14.5	
	Item No.	NTF45 070	NTF45 085	NTF45 100	NTF45 115	NTF45 130	NTF45 145
 D Ø 5.0	 7.0	 8.5	 10.0	 11.5	 13.0	 14.5	
	Item No.	NTF50 070	NTF50 085	NTF50 100	NTF50 115	NTF50 130	NTF50 145

FIXTURE

CAD/CAM SOLUTION

URIS SOLUTION FOR CEREC®

PROSTHETICS

COMPONENTS

INSTRUMENTS

1. URIS Scan Body
2. URIS Base
3. Digital Lab Analog
4. Multi-Unit Digital Components
5. Surgical Guide Sleeve



URIS Base

*2 Screws Included

URIS Scan Body

Platform	Diameter				
N Narrow Hex	● D Ø 4.0 Cuff/Height Item No.	1.0/3.5 UBN4013HA	1.0/5.5 UBN4015HA	2.0/3.5 UBN4023HA	2.0/5.5 UBN4025HA
	● D Ø 4.0 Cuff/Height Item No.	1.0/3.5 UBN4013NA	1.0/5.5 UBN4015NA	2.0/3.5 UBN4023NA	2.0/5.5 UBN4025NA
R Regular Hex	● D Ø 4.5 Cuff/Height Item No.	1.0/3.5 UBR4313HA	1.0/5.5 UBR4315HA	2.0/3.5 UBR4323HA	2.0/5.5 UBR4325HA
	● D Ø 4.5 Cuff/Height Item No.	1.0/3.5 UBR4313NA	1.0/5.5 UBR4315NA	2.0/3.5 UBR4323NA	2.0/5.5 UBR4325NA



UNSB45105H



URSB50103H

Digital Lab Analog



UDAG35D



UDAG40D

Multi-Unit Digital Components

Multi-Unit Base

*1 Screw Included



UMUB50S



UMUB50L

Multi-Unit Conversion Scan Body



UMCSB50

Multi-Unit Scan Body



UMSB50

Multi-Unit Digital Analog



UMDLA50

Surgical Guide Sleeves

N



NS38

R



NS53GR



FIXTURE

CAD/CAM SOLUTION

URIS SOLUTION FOR CEREC®

PROSTHETICS

COMPONENTS

INSTRUMENTS

1. URIS Scan Post for CEREC®
2. URIS Base for CEREC®
3. Alpha and Beta Kit



URIS Base and Scan Post for CEREC®

*URIS Base: 2 Screws Included
*Scanpost: 1 Screw and 1 Cap Included

Platform		Diameter	Gingival Height			
			GH1	GH2	GH3	GH4
Base	N Narrow Hex	D Ø 4.0 Item No.	UBN4008HS	UBN4018HS	UBN4028HS	UBN4038HS
	R Regular Hex	D Ø 4.5 Item No.	UBR4508HS	UBR4518HS	UBR4528HS	UBR4538HS
Scan Post	N Narrow Hex	D Ø 4.0 Item No.	USPN4008	USPN4018	USPN4028	USPN4038
	R Regular Hex	D Ø 4.5 Item No.	USPR4508	USPR4518	USPR4528	USPR4538

Alpha & Beta Kit

Alpha Kit

N UON-ALPHA

R UOR-ALPHA


TruScan Cap
x4


Screw
x4


GH1
x1


GH2
x1


GH3
x1


GH4
x1

Beta Kit

N UON-BETA

R UOR-BETA


TruScan Cap
x12


Screw
x12


GH1
x3


GH2
x3


GH3
x3


GH4
x3

FIXTURE

CAD/CAM SOLUTION

URIS SOLUTION FOR CEREC®

PROSTHETICS

COMPONENTS


INSTRUMENTS

1. Healing Abutment
2. Temporary Abutment
3. Direct Abutment
4. Basic Abutment
5. Milling Abutment
6. Angled Abutment
7. Multi-Unit Abutment
8. T-L Abutment
9. Ball Abutment



Healing Abutment

Hex driver: 1.27
Torque: 5-10N·cm

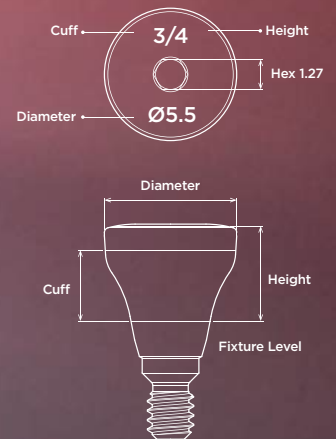

Healing Abutment



Scale 150%

Unit : mm

Platform	Diameter						
N Narrow	 D Ø 4.0						
	Cuff/Height Item No.	1.0/2.0 UHAN4012	2.0/3.0 UHAN4023	3.0/4.0 UHAN4034	4.0/5.0 UHAN4045	5.0/7.0 UHAN4057	2.0/2.0 UHAN4022
R Regular	 D Ø 4.5						
	Cuff/Height Item No.	1.0/2.0 UHAN4512	2.0/3.0 UHAN4523	3.0/4.0 UHAN4534	4.0/5.0 UHAN4545	5.0/7.0 UHAN4557	2.0/2.0 UHAN4522
	 D Ø 4.0						
	Cuff/Height Item No.	1.0/2.0 UHA4012	2.0/3.0 UHA4023	3.0/4.0 UHA4034	4.0/5.0 UHA4045	5.0/7.0 UHA4057	2.0/2.0 UHA4022
	 D Ø 4.5						
	Cuff/Height Item No.	1.0/2.0 UHA4512	2.0/3.0 UHA4523	3.0/4.0 UHA4534	4.0/5.0 UHA4545	5.0/7.0 UHA4557	2.0/2.0 UHA4522

Healing Abutment

Platform	Diameter						
R Regular	 D Ø 5.5	 1.0/2.0 UHA5512	 2.0/3.0 UHA5523	 3.0/4.0 UHA5534	 4.0/5.0 UHA5545	 5.0/7.0 UHA5557	 2.0/2.0 UHA5522
	 D Ø 6.5	 1.0/2.0 UHA6512	 2.0/3.0 UHA6523	 3.0/4.0 UHA6534	 4.0/5.0 UHA6545	 5.0/7.0 UHA6557	 2.0/2.0 UHA6522
	 D Ø 7.5	 1.0/2.0 UHA7512	 2.0/3.0 UHA7523	 3.0/5.0 UHA7535			 2.0/2.0 UHA7522

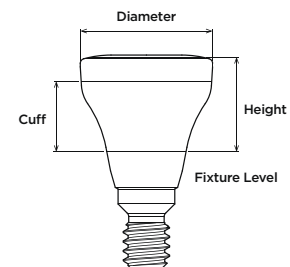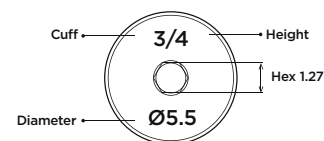

Temporary Abutment

Hex Type

Hex driver: 1.27

Torque: 20N-cm


*1 Abutment Screw Included

Non-Hex Type

Hex driver: 1.27

Torque: 20N-cm

*1 Abutment Screw Included

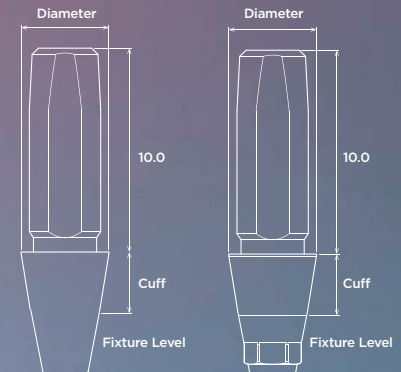



Temporary Abutment

Hex Type

Non-Hex Type


*1 Abutment Screw Included

Platform	Diameter		
Scale 150% Unit : mm N Narrow Hex	● D Ø 3.7 Cuff Item No.	● 1.0 UTAN3710H	● 3.0 UTAN3730H
	● D Ø 3.7 Cuff Item No.	● 1.0 UTAN3710N	● 3.0 UTAN3730N
R Regular Hex	● D Ø 4.3 Cuff Item No.	● 1.0 UTA4310H	● 3.0 UTA4330H
	● D Ø 4.3 Cuff Item No.	● 1.0 UTA4310N	● 3.0 UTA4330N
R Regular Non-Hex	● D Ø 4.3 Cuff Item No.		

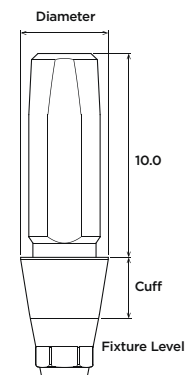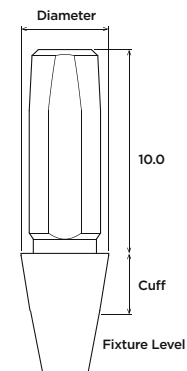

Direct Abutment

Hex driver: 1.27
Torque: Narrow [20N·cm] / Regular [30N·cm]

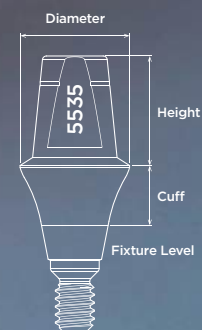

Direct Abutment

Platform

Diameter

Scale 150%

Unit : mm

N


Narrow

D Ø 4.0

Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDAN4014	UDAN4024	UDAN4034	UDAN4044	UDAN4054	UDAN4064
5.5		UDAN4015	UDAN4025	UDAN4035	UDAN4045	UDAN4055	UDAN4065
7.0		UDAN4017	UDAN4027	UDAN4037	UDAN4047	UDAN4057	UDAN4067

D Ø 4.5

Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDAN4514	UDAN4524	UDAN4534	UDAN4544	UDAN4554	UDAN4564
5.5		UDAN4515	UDAN4525	UDAN4535	UDAN4545	UDAN4555	UDAN4565
7.0		UDAN4517	UDAN4527	UDAN4537	UDAN4547	UDAN4557	UDAN4567

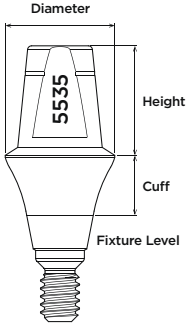

Direct Abutment

Scale 150%

Unit : mm

R

Regular

D Ø 4.5



Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDA4514	UDA4524	UDA4534	UDA4544	UDA4554	UDA4564
5.5		UDA4515	UDA4525	UDA4535	UDA4545	UDA4555	UDA4565
7.0		UDA4517	UDA4527	UDA4537	UDA4547	UDA4557	UDA4567

D Ø 5.5



Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDA5514	UDA5524	UDA5534	UDA5544	UDA5554	UDA5564
5.5		UDA5515	UDA5525	UDA5535	UDA5545	UDA5555	UDA5565
7.0		UDA5517	UDA5527	UDA5537	UDA5547	UDA5557	UDA5567

Direct Abutment

Scale 150%

Unit : mm

Platform

Diameter

D Ø 6.5

6514

6524


6534

6544

6554

6564

Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDA6514	UDA6524	UDA6534	UDA6544	UDA6554	UDA6564
5.5		UDA6515	UDA6525	UDA6535	UDA6545	UDA6555	UDA6565
7.0		UDA6517	UDA6527	UDA6537	UDA6547	UDA6557	UDA6567

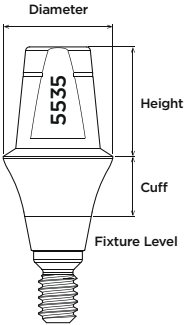

Basic Abutment

Hex Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]



*1 Abutment Screw Included

Non-Hex Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]

*1 Abutment Screw Included

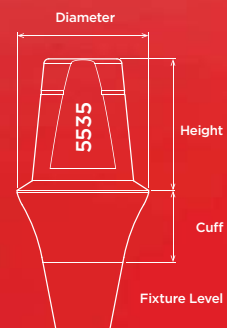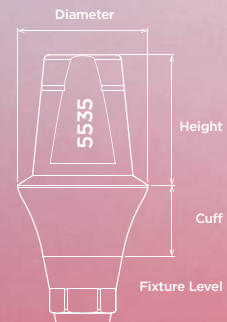

Basic Abutment

Hex Type


*1 Abutment Screw Included

Scale 150%
Unit : mm

N


Narrow
Hex

D Ø 4.0

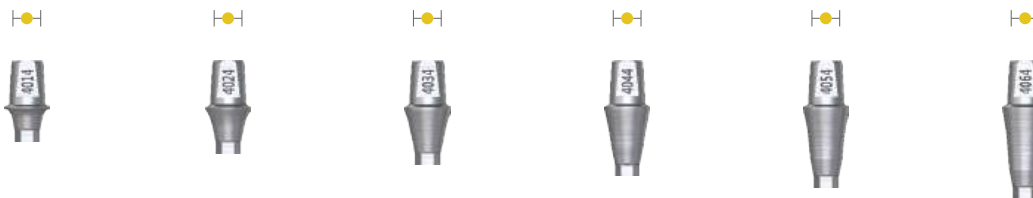


Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDAN4014H	UDAN4024H	UDAN4034H	UDAN4044H	UDAN4054H	UDAN4064H
5.5		UDAN4015H	UDAN4025H	UDAN4035H	UDAN4045H	UDAN4055H	UDAN4065H
7.0		UDAN4017H	UDAN4027H	UDAN4037H	UDAN4047H	UDAN4057H	UDAN4067H

D Ø 4.5

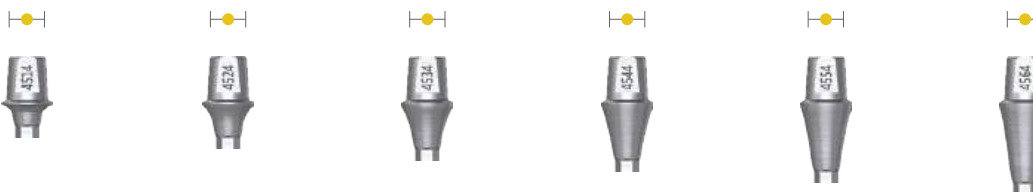

Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDAN4514H	UDAN4524H	UDAN4534H	UDAN4544H	UDAN4554H	UDAN4564H
5.5		UDAN4515H	UDAN4525H	UDAN4535H	UDAN4545H	UDAN4555H	UDAN4565H
7.0		UDAN4517H	UDAN4527H	UDAN4537H	UDAN4547H	UDAN4557H	UDAN4567H

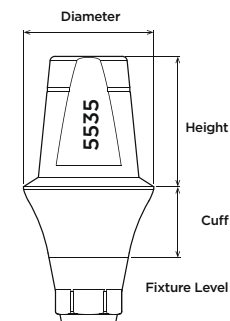

Hex Type

***1 Abutment Screw Included**

Unit: mm



Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDA4514H	UDA4524H	UDA4534H	UDA4544H	UDA4554H	UDA4564H
5.5		UDA4515H	UDA4525H	UDA4535H	UDA4545H	UDA4555H	UDA4565H
7.0		UDA4517H	UDA4527H	UDA4537H	UDA4547H	UDA4557H	UDA4567H

D Ø 5.5



Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDA5514H	UDA5524H	UDA5534H	UDA5544H	UDA5554H	UDA5564H
5.5		UDA5515H	UDA5525H	UDA5535H	UDA5545H	UDA5555H	UDA5565H
7.0		UDA5517H	UDA5527H	UDA5537H	UDA5547H	UDA5557H	UDA5567H

Basic Abutment


Hex Type

*1 Abutment Screw Included

Scale 150%

Unit : mm

Platform	Diameter								
R Regular Hex	D Ø 6.5								
		Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
		4.0		UDA6514H	UDA6524H	UDA6534H	UDA6544H	UDA6554H	UDA6564H
		5.5		UDA6515H	UDA6525H	UDA6535H	UDA6545H	UDA6555H	UDA6565H
7.0		UDA6517H	UDA6527H	UDA6537H	UDA6547H	UDA6557H	UDA6567H		

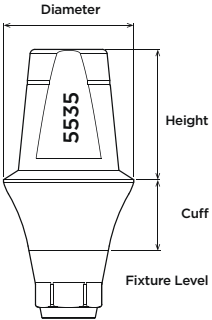

Basic Abutment Non-Hex Type

*1 Abutment Screw Included

Scale 150%


Unit : mm

N

Narrow
Non-Hex

Diameter								
 D Ø 4.0								
								
	Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
	4.0		UDAN4014N	UDAN4024N	UDAN4034N	UDAN4044N	UDAN4054N	UDAN4064N
	5.5		UDAN4015N	UDAN4025N	UDAN4035N	UDAN4045N	UDAN4055N	UDAN4065N
	7.0		UDAN4017N	UDAN4027N	UDAN4037N	UDAN4047N	UDAN4057N	UDAN4067N

 D Ø 4.5								
								
	Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
	4.0		UDAN4514N	UDAN4524N	UDAN4534N	UDAN4544N	UDAN4554N	UDAN4564N
	5.5		UDAN4515N	UDAN4525N	UDAN4535N	UDAN4545N	UDAN4555N	UDAN4565N
	7.0		UDAN4517N	UDAN4527N	UDAN4537N	UDAN4547N	UDAN4557N	UDAN4567N

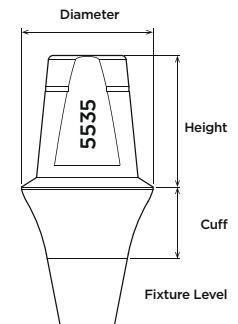

Basic Abutment

Non-Hex Type

*1 Abutment Screw Included

Platform

Diameter

Scale 150%

Unit : mm

R

Regular
Non-Hex

D Ø 4.5

4514

4524

4534

4544

4554

4564

Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDA4514N	UDA4524N	UDA4534N	UDA4544N	UDA4554N	UDA4564N
5.5		UDA4515N	UDA4525N	UDA4535N	UDA4545N	UDA4555N	UDA4565N
7.0		UDA4517N	UDA4527N	UDA4537N	UDA4547N	UDA4557N	UDA4567N

D Ø 5.5

5514

5524

5534

5544

5554

5564

Height	Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0		UDA5514N	UDA5524N	UDA5534N	UDA5544N	UDA5554N	UDA5564N
5.5		UDA5515N	UDA5525N	UDA5535N	UDA5545N	UDA5555N	UDA5565N
7.0		UDA5517N	UDA5527N	UDA5537N	UDA5547N	UDA5557N	UDA5567N


Basic Abutment

Non-Hex Type

*1 Abutment Screw Included

Scale 150%
Unit : mm

Platform	Diameter							
R Regular Non-Hex	D Ø 6.5	6514	6524	6534	6544	6554	6564	
		Cuff	1.0	2.0	3.0	4.0	5.0	6.0
4.0			UDA6514N	UDA6524N	UDA6534N	UDA6544N	UDA6554N	UDA6564N
5.5			UDA6515N	UDA6525N	UDA6535N	UDA6545N	UDA6555N	UDA6565N
7.0			UDA6517N	UDA6527N	UDA6537N	UDA6547N	UDA6557N	UDA6567N

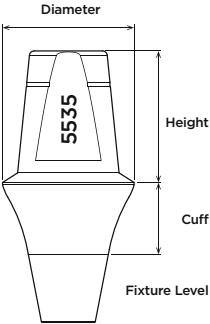

Milling Abutment

Hex Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]



*1 Abutment Screw Included

Non-Hex Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]

*1 Abutment Screw Included

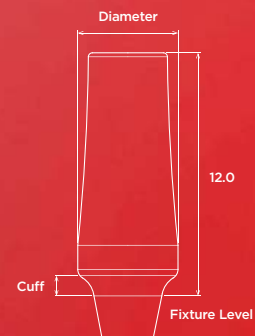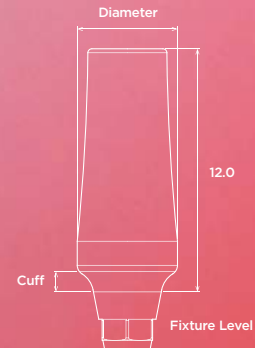

Milling Abutment

Hex Type

Non-Hex Type

*1 Abutment Screw Included

Scale 150%

Unit : mm

Platform	Diameter			Diameter		
 Narrow Hex	 D Ø 4.0 Cuff Item No.	  1.0 UMIN4010H	  3.0 UMIN4030H			
 Regular Hex	 D Ø 4.0 Cuff Item No.	  1.0 UMI4010H	  3.0 UMI4030H	 D Ø 5.0 Cuff Item No.	  1.0 UMI5010H	  3.0 UMI5030H
 Narrow Non-Hex	 D Ø 4.0 Cuff Item No.	  1.0 UMIN4010N	  3.0 UMIN4030N			
 Regular Non-Hex	 D Ø 4.0 Cuff Item No.	  1.0 UMI4010N	  3.0 UMI4030N	 D Ø 5.0 Cuff Item No.	  1.0 UMI5010N	  3.0 UMI5030N


Milling Abutment

Hex Type

Non-Hex Type

*1 Abutment Screw Included

Platform	Diameter			Diameter	Cuff	
Scale 150% Unit : mm Regular Hex	 D Ø 6.0 Cuff Item No.	 1.0 UMI6010H	 3.0 UMI6030H	 D Ø 7.0 Cuff Item No.	 1.0 UMI7010H	 3.0 UMI7030H
	 D Ø 6.0 Cuff Item No.	 1.0 UMI6010N	 3.0 UMI6030N	 D Ø 7.0 Cuff Item No.	 1.0 UMI7010N	 3.0 UMI7030N

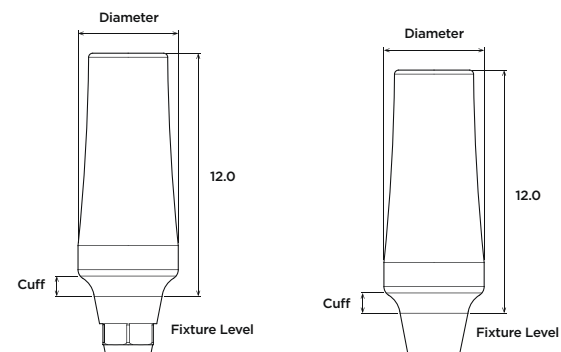

Angled Abutment 17°

Hex Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]



*1 Abutment Screw Included

Non-Hex Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]

*1 Abutment Screw Included

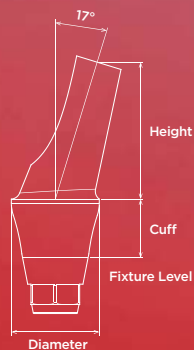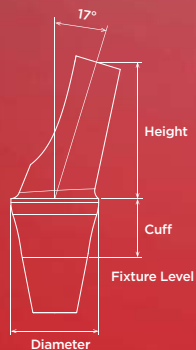

Angled Abutment 17° Hex Type

*1 Abutment Screw Included

Platform	Diameter					
Scale 150% Unit : mm Narrow Hex	D Ø 4.0 Cuff Item No.					
						
						
		2.0 U17AAN4027H	3.0 U17AAN4037H	4.0 U17AAN4047H	5.0 U17AAN4057H	
Regular Hex	D Ø 4.5 Cuff Item No.					
						
						
		2.0 U17AA4527H	3.0 U17AA4537H	4.0 U17AA4547H	5.0 U17AA4557H	
	D Ø 5.5 Cuff Item No.					
						
						
		2.0 U17AA5527H	3.0 U17AA5537H	4.0 U17AA5547H	5.0 U17AA5557H	

Angled Abutment 17°


Non-Hex Type

*1 Abutment Screw Included

Scale 150%

Unit : mm

Platform	Diameter				
N Narrow Hex	D Ø 4.0				
					
			2.0	3.0	4.0
		Cuff Item No.	U17AAN4027N	U17AAN4037N	U17AAN4047N
R Regular Non-Hex	D Ø 4.5				
					
			2.0	3.0	4.0
		Cuff Item No.	U17AA4527N	U17AA4537N	U17AA4547N
R Regular Non-Hex	D Ø 5.5				
					
			2.0	3.0	4.0
		Cuff Item No.	U17AA5527N	U17AA5537N	U17AA5547N

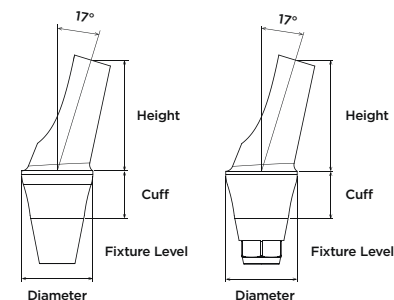

Multi-Unit Abutment

Straight Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]



*1 Holder Included

Angled Type

Hex driver: 1.27

Torque: Narrow [20N·cm] / Regular [30N·cm]

*1 Abutment Screw and 1 Driver Included

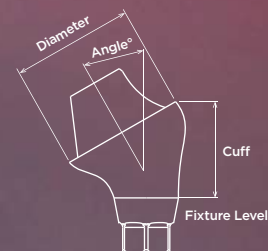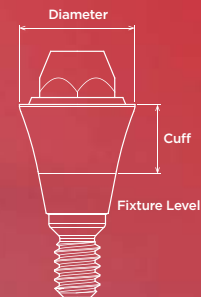

Abutment Screw

Multi-Unit Abutment


Straight Type

*1 Holder Included

Scale 150%

Unit : mm


Platform	Diameter						
N Narrow							
	D Ø 5.0						
	Cuff Item No.	1.0 UMAN5010	2.0 UMAN5020	3.0 UMAN5030	4.0 UMAN5040	5.0 UMAN5050	6.0 UMAN5060
R Regular							
	D Ø 5.0						
	Cuff Item No.	1.0 UMA5010	2.0 UMA5020	3.0 UMA5030	4.0 UMA5040	5.0 UMA5050	6.0 UMA5060

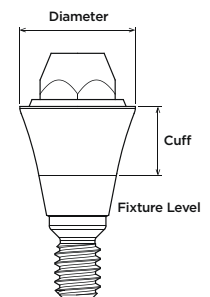

Multi-Unit Abutment Angled Type

*1 Abutment Screw and 1 Driver Included

Platform	Diameter						
Scale 150% Unit : mm Narrow Hex	 D Ø 5.0	 	 	 	 	 	 
	Cuff/Angle Item No.	3.0/17° U17MUAN5030H	4.0/17° U17MUAN5040H	5.0/17° U17MUAN5050H	4.0/29.5° U30MUAN5040H	5.0/29.5° U30MUAN5050H	6.0/29.5° U30MUAN5060H
Regular Hex	 D Ø 5.0	 	 	 	 	 	 
	Cuff/Angle Item No.	3.0/17° U17MUAR5030H	4.0/17° U17MUAR5040H	5.0/17° U17MUAR5050H	4.0/29.5° U30MUAR5040H	5.0/29.5° U30MUAR5050H	6.0/29.5° U30MUAR5060H

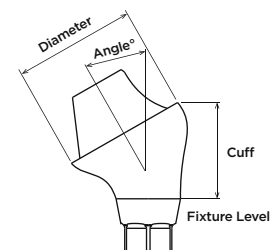



T-L Analog
Pg.75



T-L Impression Coping
Pg.75



T-L Abutments

T-L Abutment

Scale 150%

Unit : mm

Platform	Diameter						
N Narrow	 D Ø 3.8	  1.0 UTLAN3810	  2.0 UTLAN3820	  3.0 UTLAN3830	  4.0 UTLAN3840	  5.0 UTLAN3850	  6.0 UTLAN3860
	Cuff Item No.						
R Regular	 D Ø 3.8	  1.0 UTLAR3810	  2.0 UTLAR3820	  3.0 UTLAR3830	  4.0 UTLAR3840	  5.0 UTLAR3850	  6.0 UTLAR3860
	Cuff Item No.						

Ball Abutment

Straight Type

For Oververdenture Implant

Torque: Narrow [20N·cm] / Regular [30N·cm]

*1 Holder Included

O-Ring


Pg.75

Ball Abutment Retainer

Pg.75

Ball Abutment Retainer Cap


Pg.75

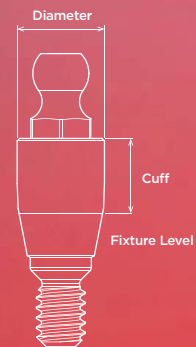

Ball Abutment


*1 Holder Included

Platform	Diameter						
Scale 150% Unit : mm N Narrow	● D Ø 3.5 Cuff Item No.						
							
		1.0 UBAN3510	2.0 UBAN3520	3.0 UBAN3530	4.0 UBAN3540	5.0 UBAN3550	6.0 UBAN3560
							
							
		1.0 UBA3510	2.0 UBA3520	3.0 UBA3530	4.0 UBA3540	5.0 UBA3550	6.0 UBA3560

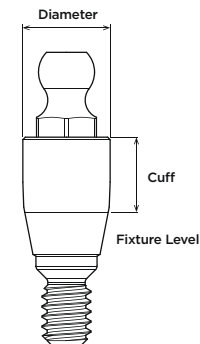

URIS Digital Solution

Powered by **TruAbutment**

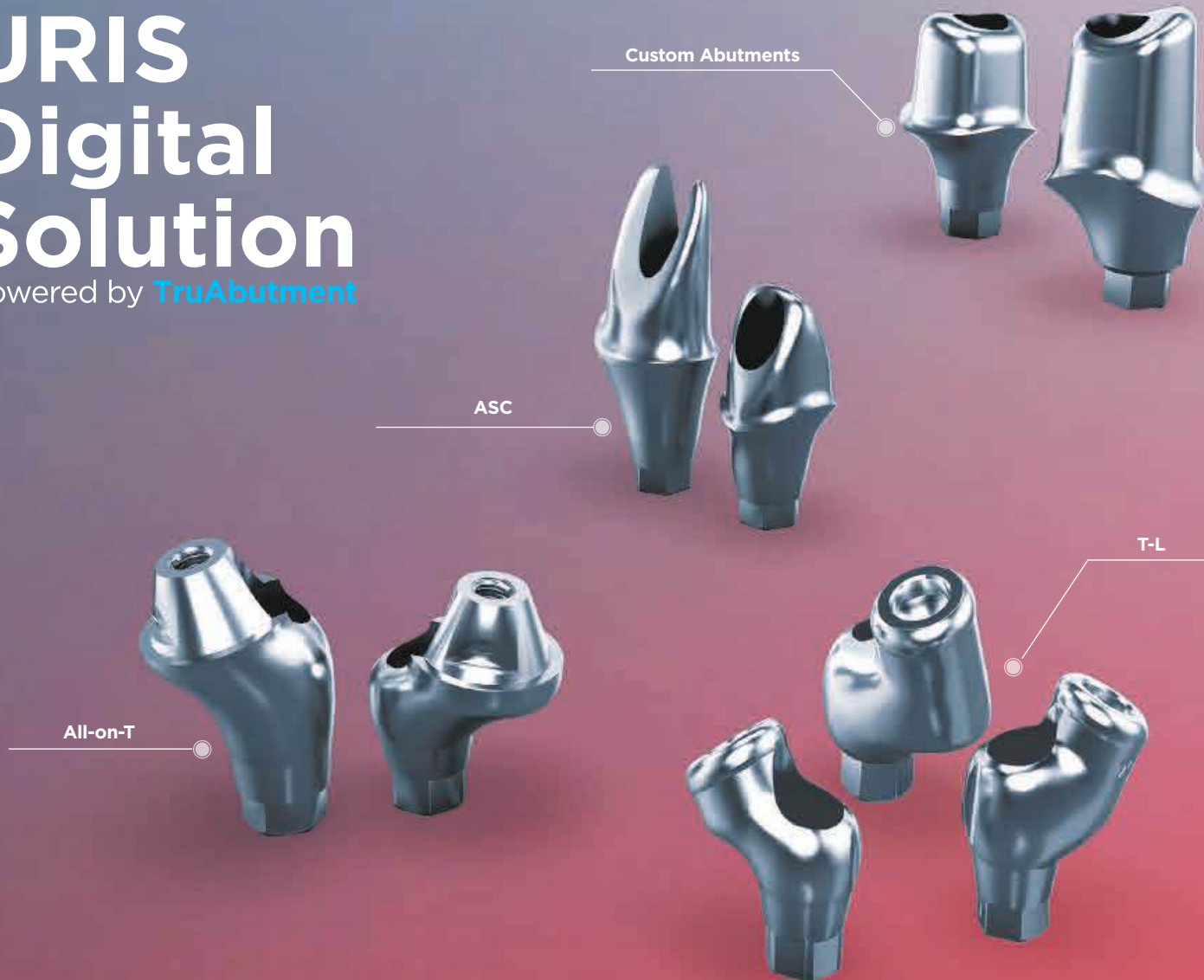

FIXTURE

CAD/CAM SOLUTION


URIS SOLUTION FOR CEREC®

PROSTHETICS

COMPONENTS





INSTRUMENTS

1. Impression Components
2. Protective Cap
3. Impression Cap
4. Abutment Components
5. Multi-Unit Components
6. T-L Components
7. Ball Components
8. Fixture Analog
9. Direct Abutment Analog
10. Ball Abutment Analog
11. Abutment Screws

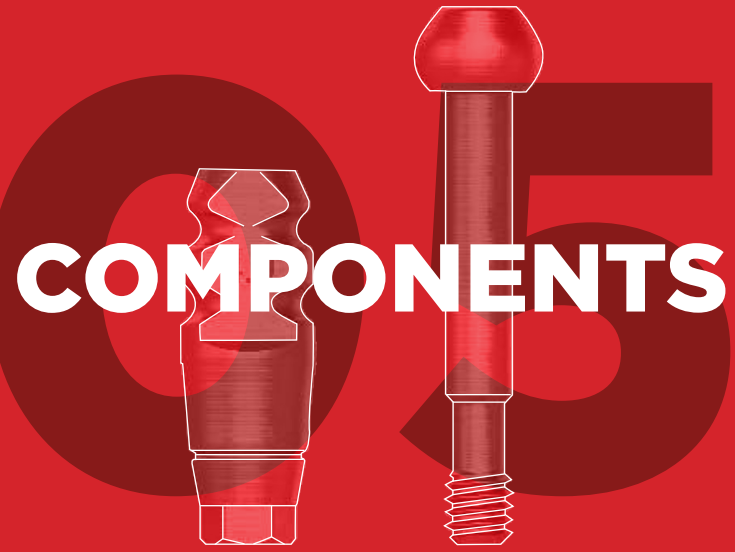

Impression Components





Transfer Type

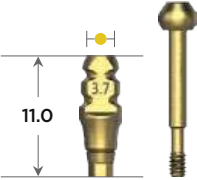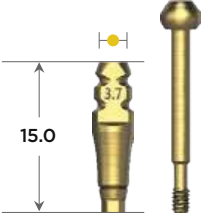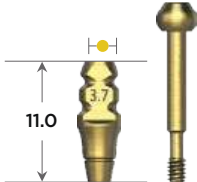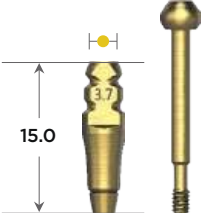Platform	Diameter		
Scale 150% Unit : mm N Narrow Hex	D Ø 3.7	 UIPTN3711H	 UIPTN3714H
	D Ø 3.7	 UIPTN3711N	 UIPTN3714N

Impression Components Transfer Type





Scale 150%

Unit : mm

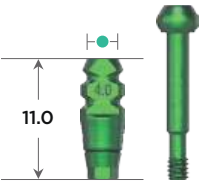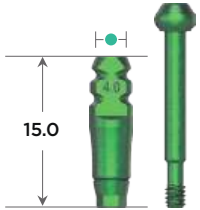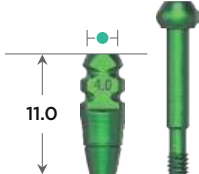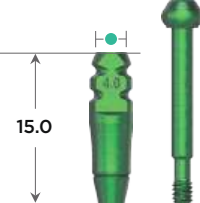Platform	Diameter		
R Regular Hex			
	D Ø 4.0	UIPT4011H	UIPT4014H
	D Ø 4.5	UIPT4511H	UIPT4514H
	D Ø 5.5	UIPT5511H	UIPT5514H
R Regular Non-Hex			
	D Ø 4.0	UIPT4011N	UIPT4014N
	D Ø 4.5	UIPT4511N	UIPT4514N
	D Ø 5.5	UIPT5511N	UIPT5514N

Impression Components





Pick Up Type

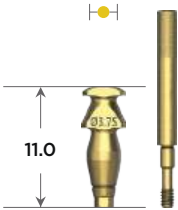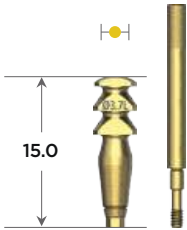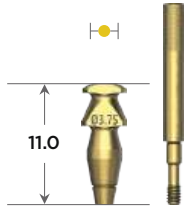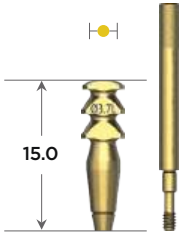Platform	Diameter		
Scale 150% Unit : mm	Narrow Hex	D Ø 3.7	
			
	Narrow Non-Hex	D Ø 3.7	
			

Impression Components Pick Up Type

Scale 150%

Unit : mm

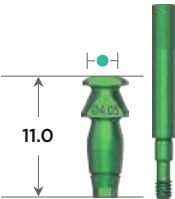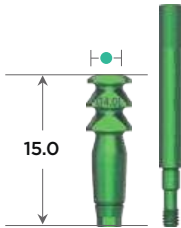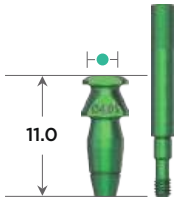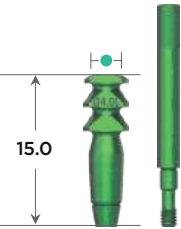Platform	Diameter		
R Regular Hex			
	D Ø 4.0	UIPP4011H	UIPP4015H
	D Ø 4.5	UIPP4511H	UIPP4515H
	D Ø 5.5	UIPP5511H	UIPP5515H
R Regular Non-Hex			
	D Ø 4.0	UIPP4011N	UIPP4015N
	D Ø 4.5	UIPP4511N	UIPP4515N
	D Ø 5.5	UIPP5511N	UIPP5515N

Protective Cap

Diameter			
D Ø 4.0			
	4.0 UDAC4004	5.5 UDAC4005	7.5 UDAC4007
	Height Item No.		
D Ø 4.5			
	4.0 UDAC4504	5.5 UDAC4505	7.5 UDAC4507
	Height Item No.		
D Ø 5.5			
	4.0 UDAC5504	5.5 UDAC5505	7.5 UDAC5507
	Height Item No.		
D Ø 6.5			
	4.0 UDAC6504	5.5 UDAC6505	7.5 UDAC6507
	Height Item No.		

Impression Cap

Ø 4.0	Ø 4.5	Ø 5.5	Ø 6.5
			
UDAIC40	UDAIC45	UDAIC55	UDAIC65

Abutment Components

For Single Case Use

Ø 4.0	Ø 4.5	Ø 5.5	Ø 6.5
			
UDABC40S	UDABC45S	UDABC55S	UDABC65S

For Bridge Case Use

Ø 4.0	Ø 4.5	Ø 5.5	Ø 6.5
			
UDABC40B	UDABC45B	UDABC55B	UDABC65B

Multi-Unit Components

Impression Post

Pick Up



UMIPP50N

Impression Post Screw

Pick Up



UMIPPS07

Impression Post

Transfer



UMIPT50

Multi-Unit Analog



UMAG50

Healing Cap



UMHC50

Multi-Unit Components

*1 Multi-Unit Cylinder Screw Included

Scale 150%

Unit : mm

Temporary Cylinder



UMTC50N

Burn Out Cylinder



UMBC50N

NEW

Multi-Unit Direct Screw



UMDW16

NEW

Multi-Unit Direct Screw Set



x 5

UMDW16(5)



x 20

UMDW16(20)

T-L Components

Impression Coping



UTLP55

Analog



UTLG38

Fixture Analog

N



UDAG35

R



UDAG40

Ball Components

*1 O-Ring Included

Retainer Cap



UBSC35

+

O-Ring



UOR4515K

Retainer



UBSC35

+

O-Ring



UOR4515K

Direct Abutment Analog

Diameter			
D Ø 4.0			
			
	Height Item No.	4.0 UDAG404	5.5 UDAG405
D Ø 4.5			
			
	Height Item No.	4.0 UDAG454	5.5 UDAG455
D Ø 5.5			
			
	Height Item No.	4.0 UDAG554	5.5 UDAG555
D Ø 6.5			
			
	Height Item No.	4.0 UDAG654	5.5 UDAG655

7.0
UDAG407

7.0
UDAG457

7.0
UDAG557

7.0
UDAG657

Ball Abutment Analog

Scale 150%

Unit : mm



UBAG35

Abutment Screws

N



UAS16H

R



UAS20H

ASC (Angulated Screw Channel) Abutment Screws

N



UON-AW

R



UOR-AW

FIXTURE

CAD/CAM SOLUTION


URIS SOLUTION FOR CEREC®

PROSTHETICS

COMPONENTS

INSTRUMENTS

1. Surgical Kit
2. Wide Surgical Kit
3. Tapered Surgical Kit
4. Surgical Plus Kit
5. Pylon Kit
6. Pylon Plus Kit
7. Tapered Pylon Kit
8. Narrow Pylon Kit
9. Pylon Anchor Kit
10. Pylon Crestal Sinus Kit
11. Prosthetic Kit
12. ASC Kit

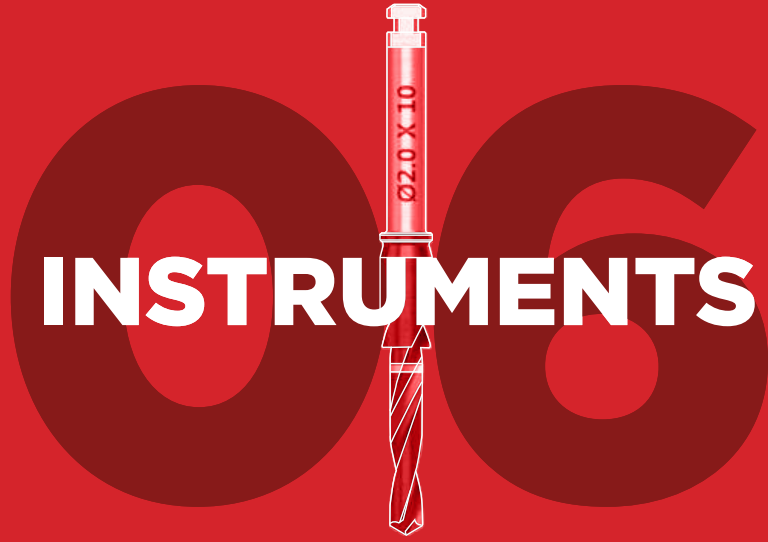


SURGICAL KIT

SKA01



Surgical Kit Item No. SKA01

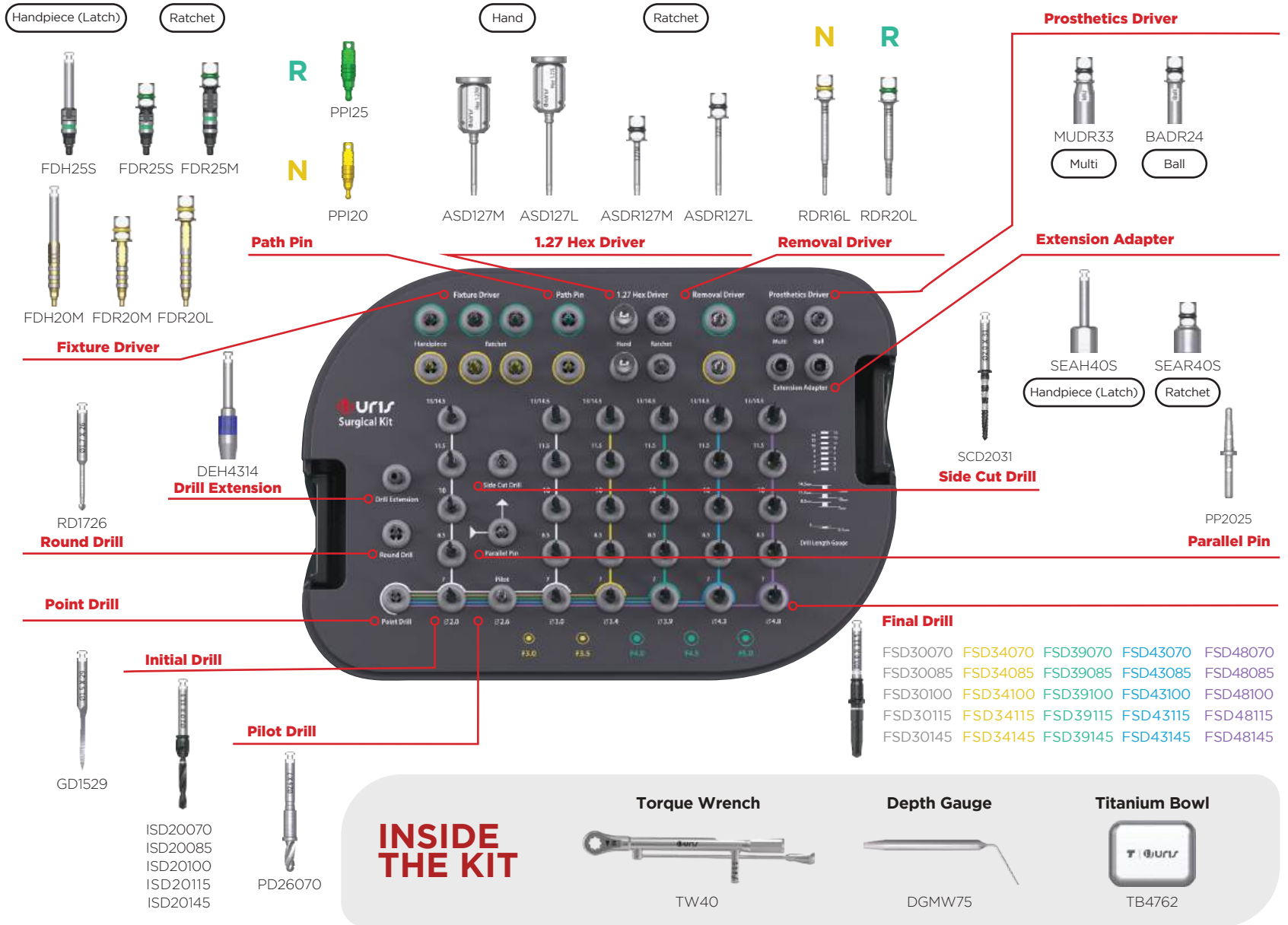

Point Drill Surgical

Drills through the cortical bone to create an ideal path for the next drills at the selected site.
800 - 1200 rpm (must use with saline water)

Scale 150%

Unit : mm



GD1529

Initial Drill Surgical

Expands the prepared site and used according to the chosen fixture's length.


800 - 1200 rpm

Unit : mm

Scale 1:1

D Ø 2.0

Item No.

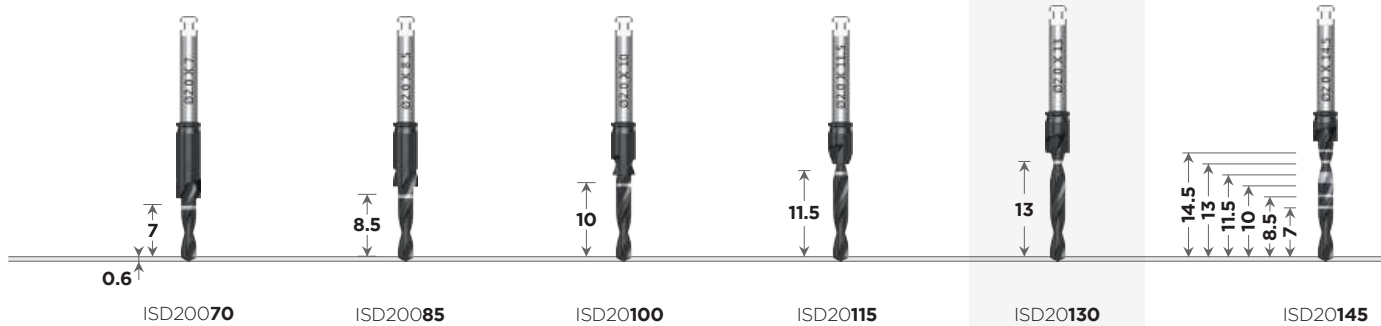

*Sold Separately

Parallel Pin Surgical

Indicates and confirms the bone preparation location.

Scale 150%

Unit : mm



PP2025

Pilot Drill Surgical

Expands cortical bone to let the final drill enter the path easily.
800 - 1200 rpm



PD26070

Final Drill Surgical


The last drill used to place the URIS fixture and is available for each fixture diameter/length and bone density. 800 - 1200 rpm

Diameter

*Sold Separately

D Ø 3.0

Item No.



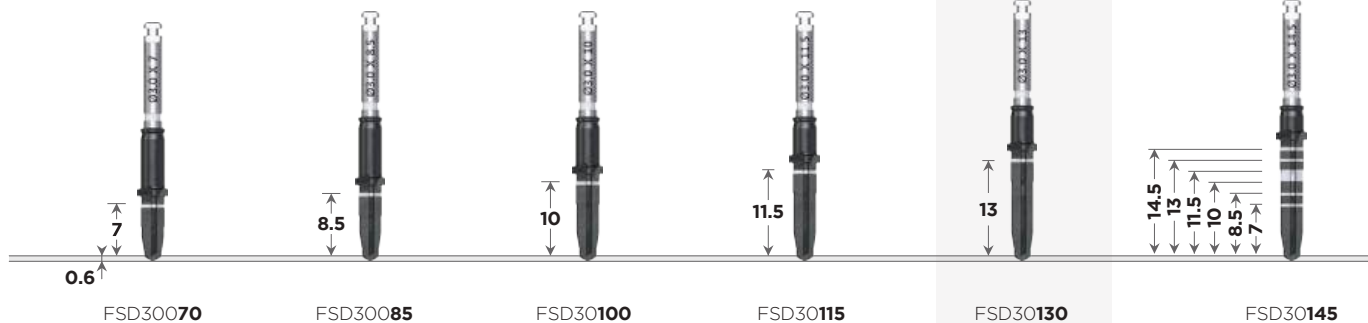


Final Drill

Surgical

The last drill used to place the URIS fixture and is available for each fixture diameter/length and bone density.

800 - 1200 rpm

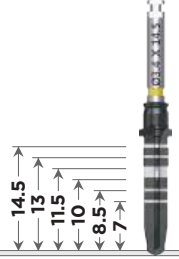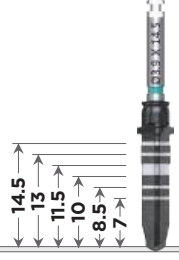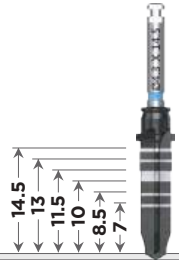Diameter							
Scale 150% Unit : mm D Ø 3.4	Item No.						
		FSD34070	FSD34085	FSD34100	FSD34115	FSD34130	FSD34145
D Ø 3.9	Item No.						
		FSD39070	FSD39085	FSD39100	FSD39115	FSD39130	FSD39145
D Ø 4.3	Item No.						
		FSD43070	FSD43085	FSD43100	FSD43115	FSD43130	FSD43145

*Sold Separately

Final Drill **Surgical** The last drill used to place the URIS fixture and is available for each fixture diameter/length and bone density. 800 - 1200 rpm

Scale 150%

Unit : mm

D Ø 4.8


Item No.

Diameter					
					
	0.7	8.5	10	11.5	13
	FSD48070	FSD48085	FSD 48100	FSD48115	FSD 48130
					
					14.5
					13
					11.5
					10
					8.5
					7
					FSD48145

*Sold Separately

Side Cut Drill **Surgical**

Drills with the side blade to modify the drilling path.
800 - 1200 rpm

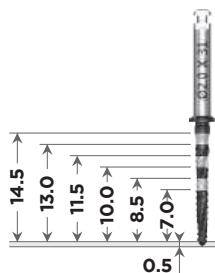

SCD2031

Path Pin **Surgical**

Confirms the path, the gingiva depth & the hex direction of the prosthetic connection after the placement.

N



PPI20

R



PPI25

Drill Extension

Extends the length of latch type tools.



DEH4314

Hex Driver 1.27

*Caution: Must be properly engaged and used at the recommended torque value.

Scale 150%

Unit : mm

Ratchet



ASDR127M



ASDR127L

Hand



ASD127M



ASD127L

Extension Adapter

Extends the instrument length.

Handpiece (Latch)



SEAH40S



SEAH40L
*Sold Separately

Ratchet



SEAR40S



SEAR40L
*Sold Separately

Ball Abutment Driver

Compatible with the ball abutment.

Ratchet



BADR24

Multi-Unit Abutment Driver

Compatible with Straight Multi-Unit abutment.

Ratchet



MUDR33

Round Drill **Surgical**

Levels out the uneven bone and to remove the remaining gingiva residue.

800 - 1200 rpm



RD1726

Removal Driver **Surgical**

Used to remove the abutment off the fixture. Needs to be inserted upright and tightened clockwise.

N

R



RDR16S

***Sold Separately**



RDR16L



RDR20S

***Sold Separately**



RDR20L

Depth Gauge

Measures the depth of a drilling hole.



DGMW75

Torque Wrench

Tightens prosthetic components with the use of compatible ratchet-type drivers.

(Use at the recommended torque value of 20 - 40 N·cm)



TW40

Titanium Bowl



TB4762

Fixture Driver

Surgical

*Caution: Must torque after fully engaging into the fixture securely to avoid any breakages.

Platform

Scale 150%

Unit : mm

N
Narrow

Handpiece (Latch)

Ratchet



FDH20S
*Sold Separately



FDH20M



FDH20L
*Sold Separately



FDR20S
*Sold Separately



FDR20M



FDR20L

R
Regular



FDH25S



FDH25M
*Sold Separately



FDH25L



FDR25S




FDR25M




FDR25L
*Sold Separately

Surgical Drilling Protocol OMNI Fixture

Implant Ø 3.0 x 11.5mm

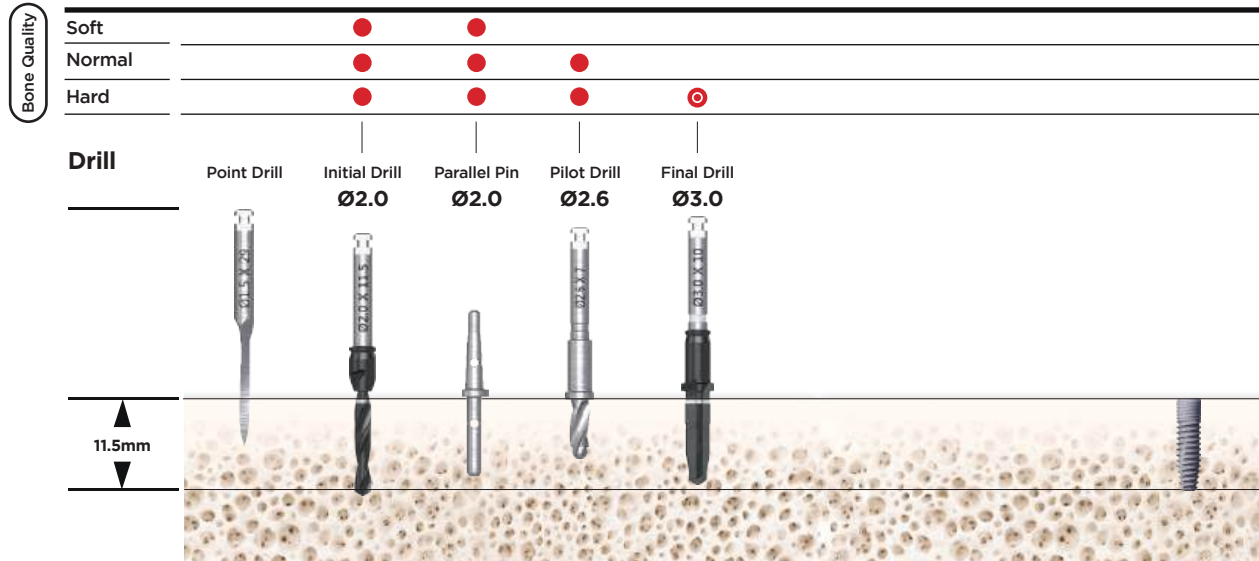


Implant Ø 3.5 x 11.5mm

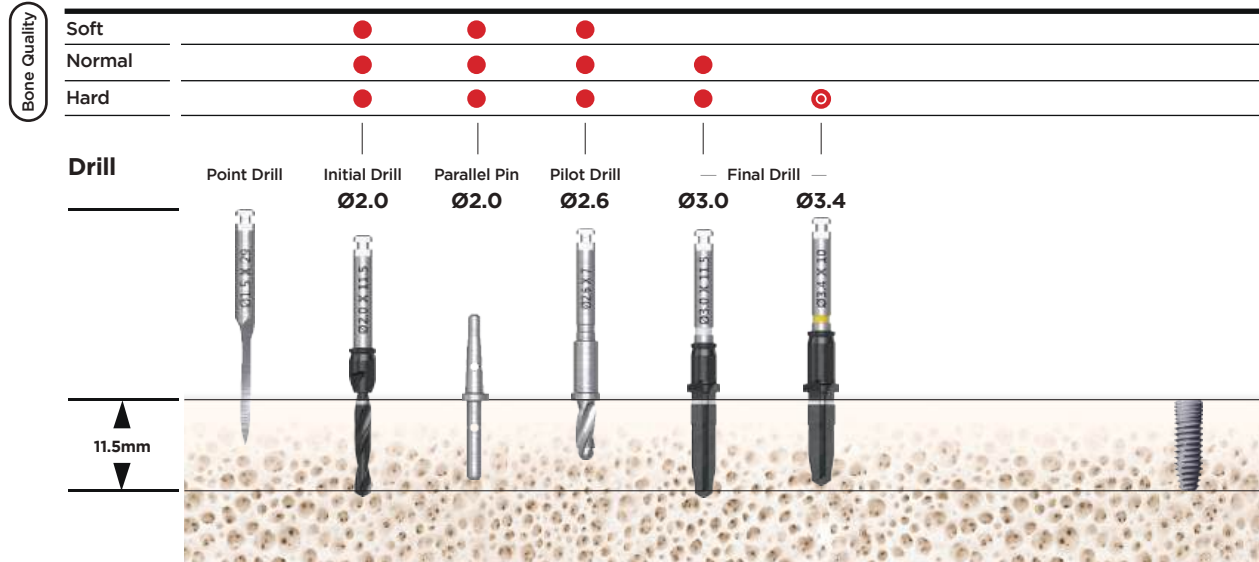


- Required
- Use at your discretion
- ⊙ Requires a drill that is 1.5mm shorter than the fixture

Surgical Drilling Protocol OMNI Fixture

Implant Ø 4.0 x 11.5mm

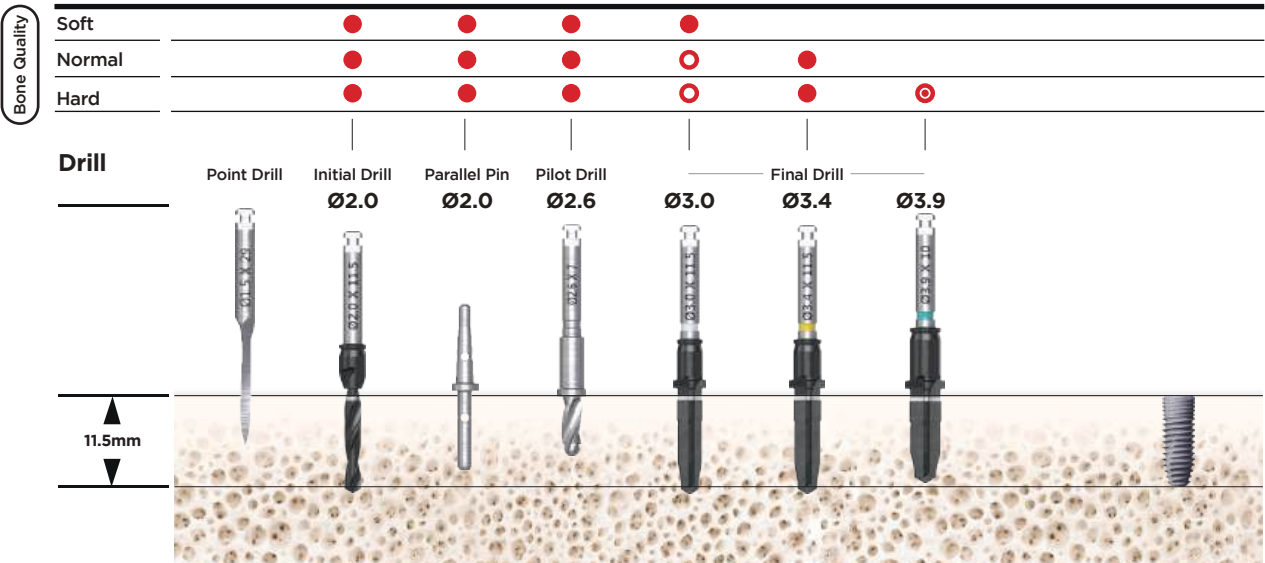


Implant Ø 4.5 x 11.5mm

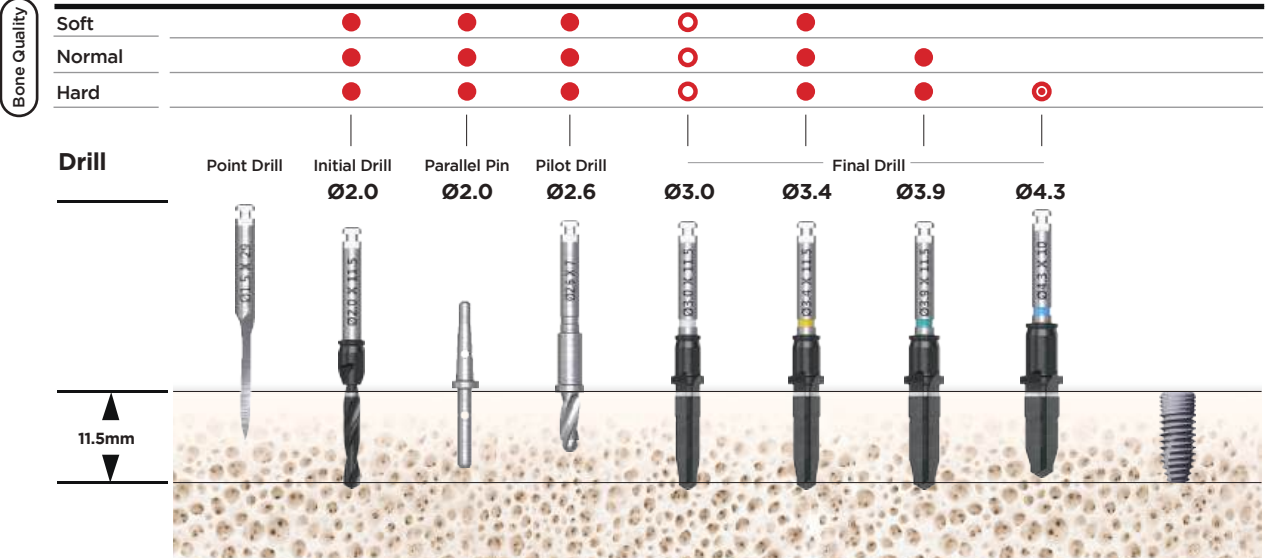

- Required
- Use at your discretion
- ⊙ Requires a drill that is 1.5mm shorter than the fixture

Implant Ø 5.0 x 11.5mm

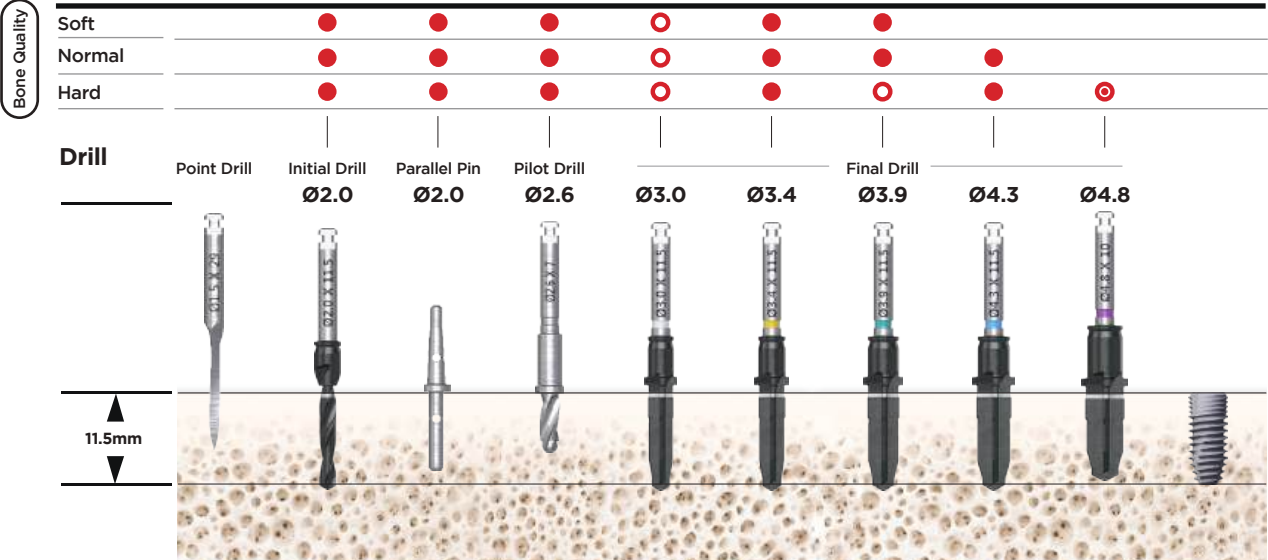

WIDE SURGICAL KIT



WSK01



Wide Surgical Kit

Item No. WSK01



Final Drill Surgical




The last drill used to place the URIS fixture and is available for each fixture diameter/length and bone density. 800 - 1200 rpm

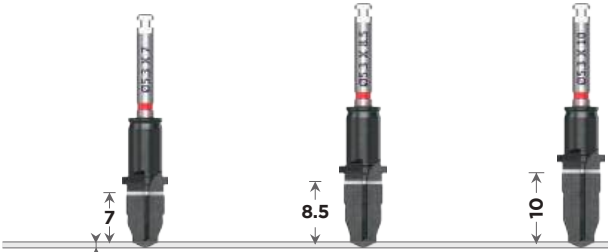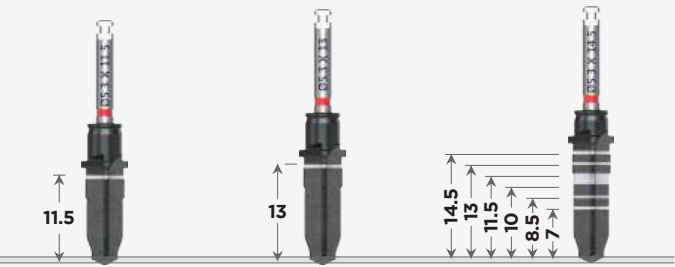Diameter		Length					
Scale 150% Unit : mm	 D Ø 5.3						
		FSD53070	FSD53085	FSD53100	*Sold Separately 		
Item No.							
Scale 150% Unit : mm	 D Ø 5.8						
		FSD58070	FSD58085	FSD58100			
Item No.							

Final Drill

Surgical


The last drill used to place the URIS fixture and is available for each fixture diameter/length and bone density. 800 - 1200 rpm

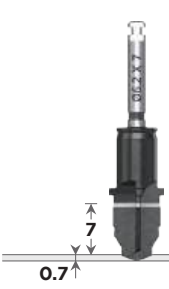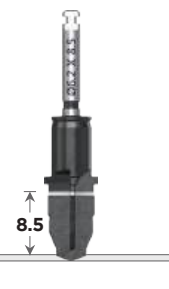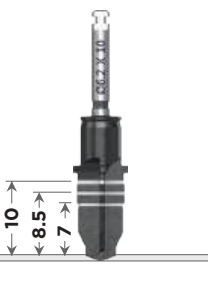Diameter	Length		
Scale 150% Unit : mm D Ø 6.2 Item No.	 FSD62070	 FSD62085	 FSD62100




Surgical Drilling Protocol OMNI Wide Fixture

Implant Ø 5.5 x 11.5mm

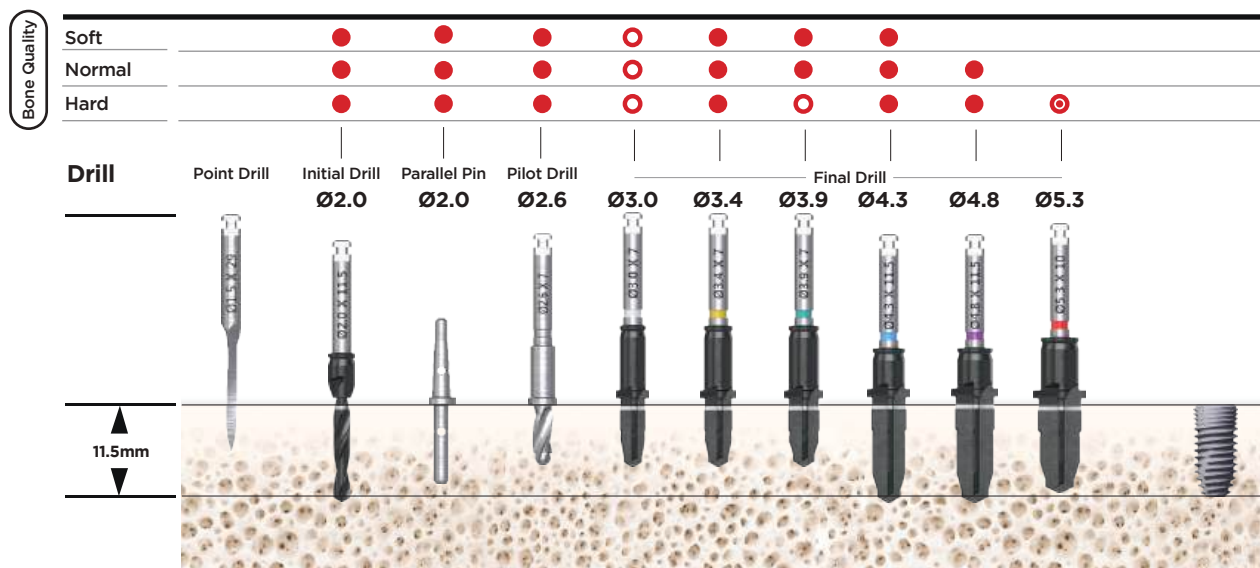


Implant Ø 6.0 x 10.0mm

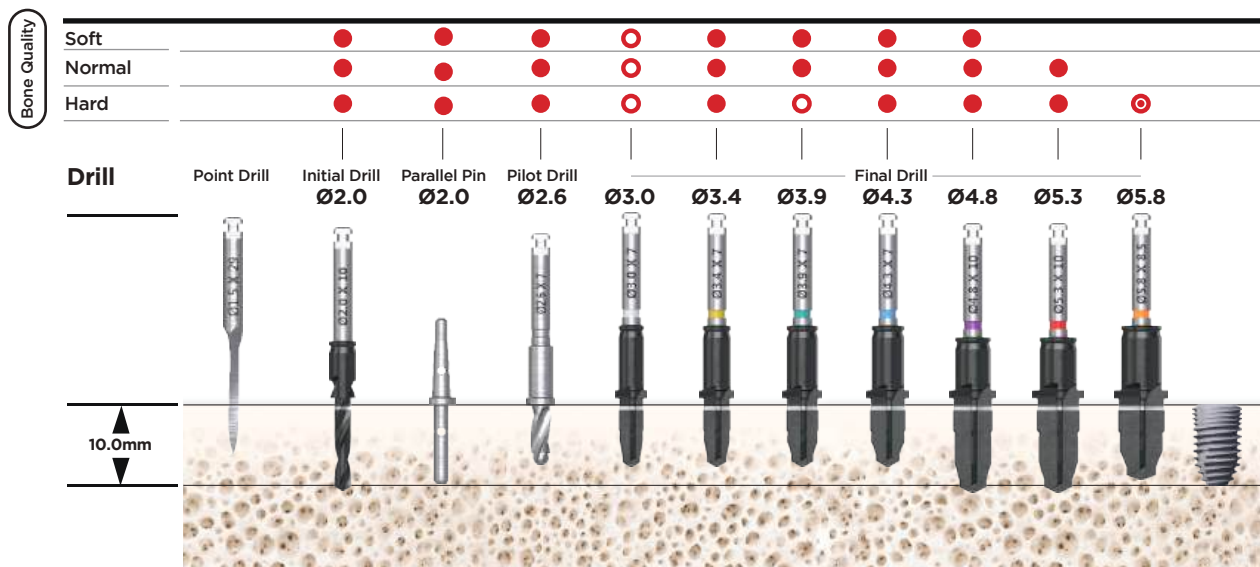

- Required
- Use at your discretion
- ⊙ Requires a drill that is 1.5mm shorter than the fixture

Implant Ø 6.5 x 10.0mm

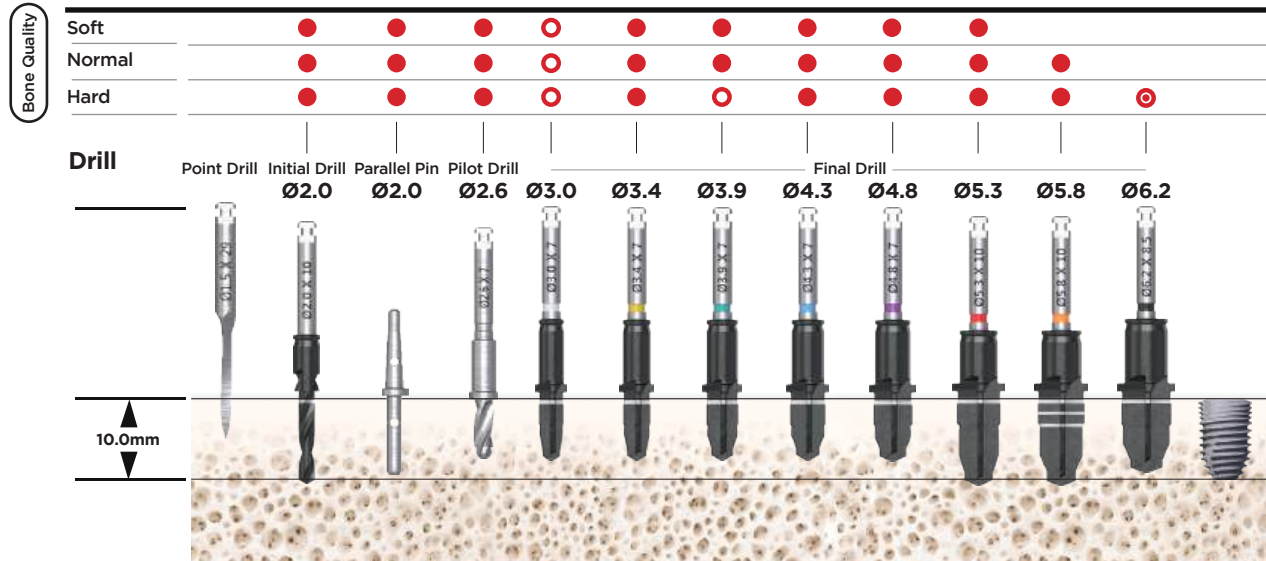


TAPERED SURGICAL KIT

TSK01



Tapered Surgical Kit

Item No. TSK01

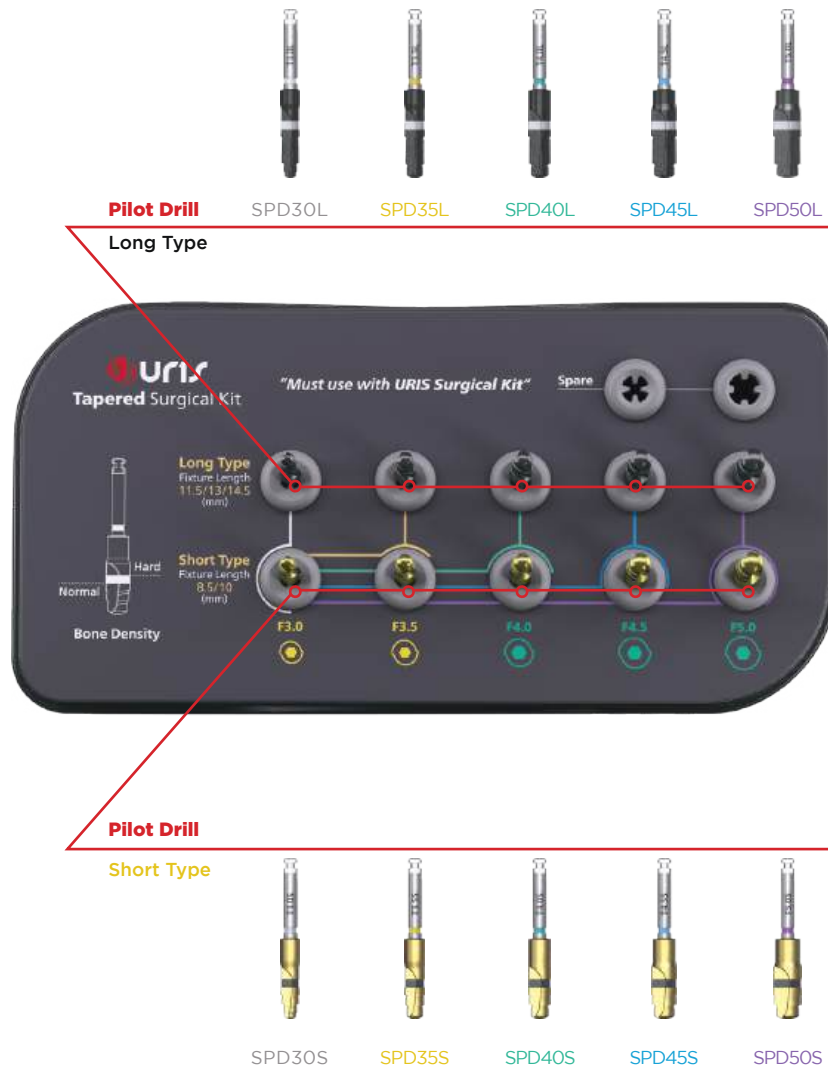

Pilot Drill

Surgical

Used as the last drill to place a URIS tapered fixture and is selected according to the tapered fixture diameter and bone density.

Low Speed: 50 - 100 rpm; within 5 seconds with high torque
High Speed: 800 - 1200 rpm (use irrigation)

Length		Diameter					
Scale 150%	Unit : mm	Long Type	Ø 3.0	Ø 3.5	Ø 4.0	Ø 4.5	Ø 5.0
							
			Item No.	SPD30L	SPD35L	SPD40L	SPD45L

Pilot Drill **Surgical**

Used as the last drill to place a URIS tapered fixture and is selected according to the tapered fixture diameter and bone density.

Low Speed: 50 - 100 rpm; within 5 seconds with high torque
High Speed: 800 - 1200 rpm (use irrigation)

Scale 150%
Unit : mm

Length		Diameter				
Short Type		Ø 3.0	Ø 3.5	Ø 4.0	Ø 4.5	Ø 5.0
						
Item No.		SPD30S	SPD35S	SPD40S	SPD45S	SPD50S




SURGICAL PLUS KIT

SPK01

Surgical Plus Kit

Item No. SPK01

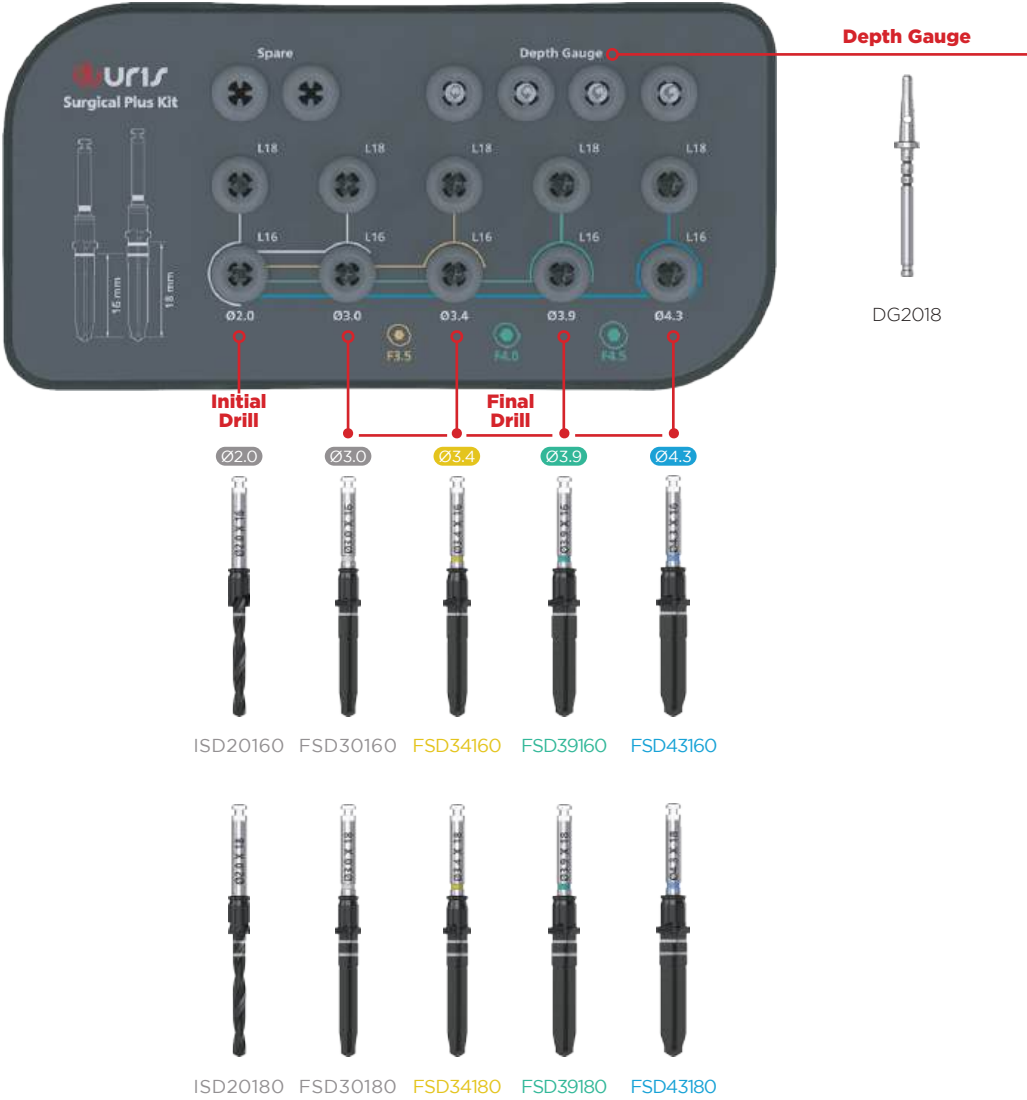

Initial Drill Surgical

Expands the prepared site and used according to the chosen fixture's length.

800 - 1200 rpm



Depth Gauge

Ø 2.0



DG2018

Final Drill

Surgical

The last drill used to place a URIS fixture and is available for each fixture diameter/length and bone density.

Diameter




Scale 150%
Unit : mm


GUIDED

PYLON KIT

NKA01

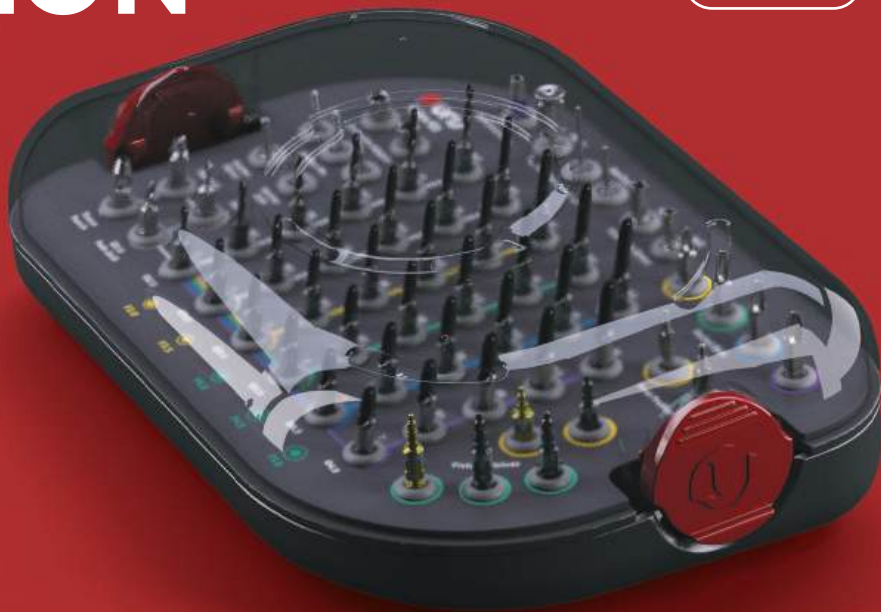

Pylon Kit Item No. NKA01

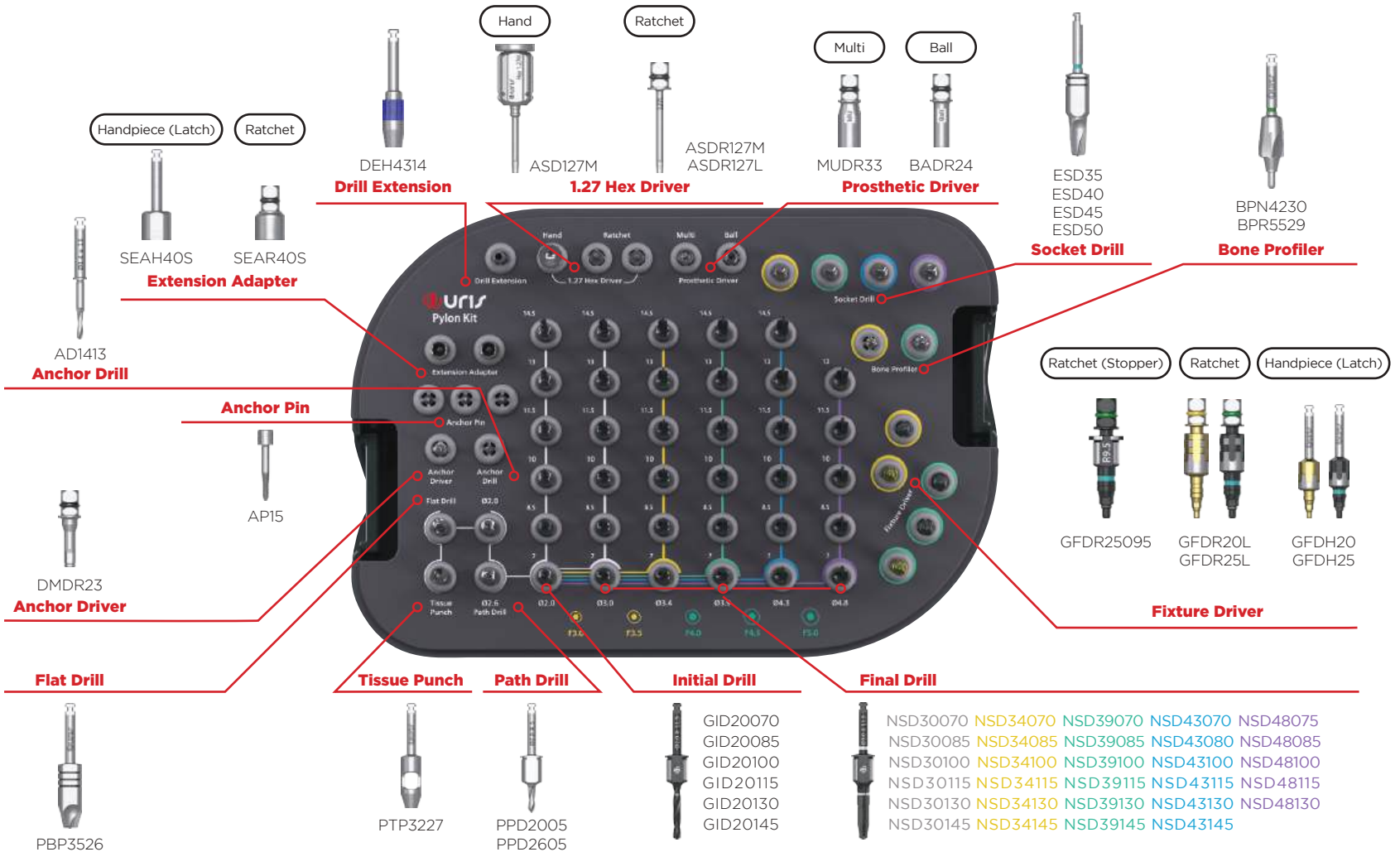

**INSIDE
THE KIT**

Torque Wrench



TW40

Titanium Bowl




TB4762

Flat Drill Pylon

Flattens out uneven bone and also removes remaining gingiva residue after tissue punching.

Scale 150%

Unit : mm

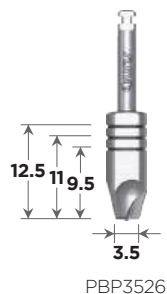


Tissue Punch Pylon

Removes soft tissue during a flap-less surgery.



Path Drill Pylon

Increases the accuracy of the drilling path.

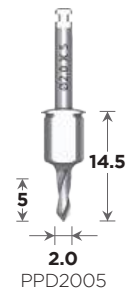

Tissue Punch/Path Drill/Initial Drill

Low speed: 50 - 100 rpm; within 5 seconds with high torque | High Speed: 800 - 1200 rpm (use irrigation)


Initial Drill Pylon

Expands the drilling path by the chosen fixture's length.

Diameter

D Ø 2.0

Item No.

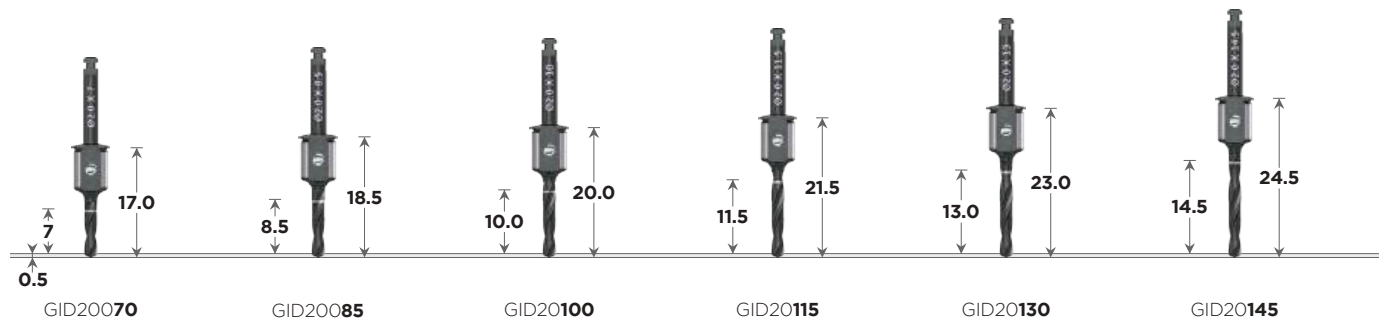

Final Drill

Pylon

Used as the last drill to place a URIS fixture and is selected according to the fixture diameter and bone density.


Low speed: 50 - 100 rpm; within 5 seconds with high torque
High speed: 800 - 1200 rpm (use irrigation)

Scale 150%

Unit : mm


D Ø 3.0

Item No.

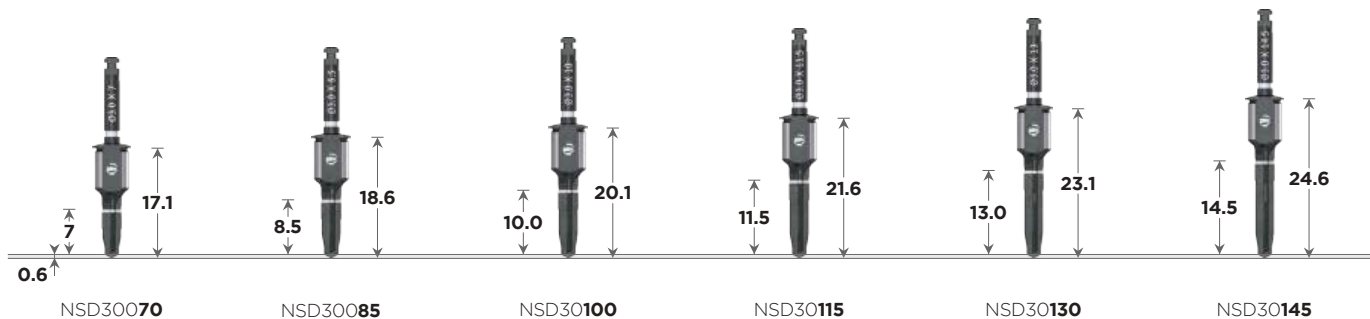


D Ø 3.4

Item No.

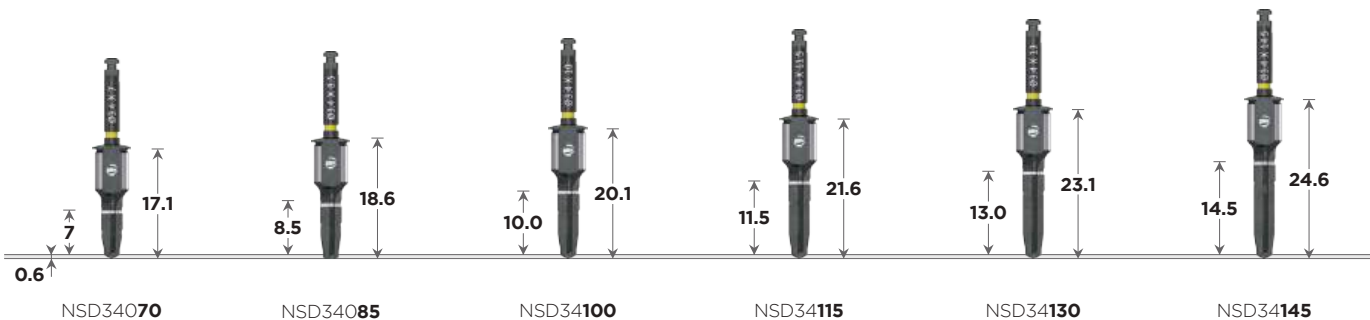

D Ø 3.9







Item No.

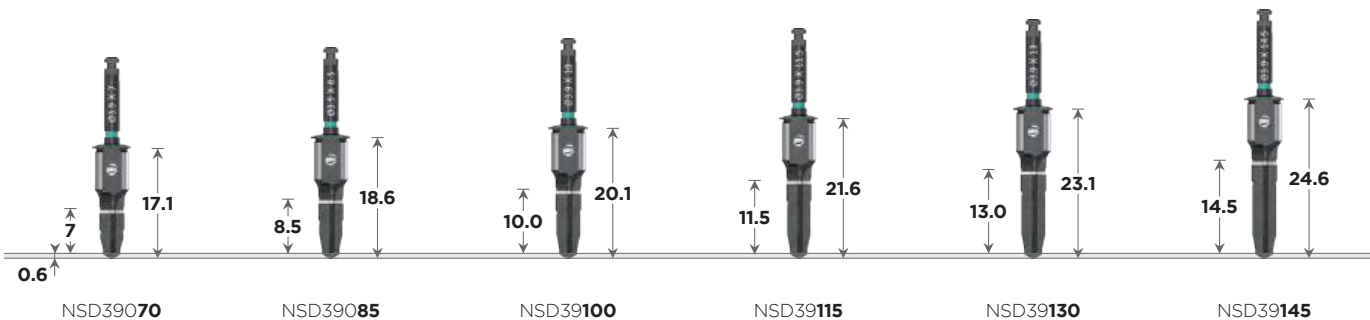

Final Drill Pylon

Used as the last drill to place a URIS fixture and is selected according to the fixture diameter and bone density.

Low speed: 50 - 100 rpm; within 5 seconds with high torque
High speed: 800 - 1200 rpm (use irrigation)

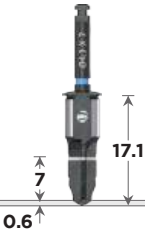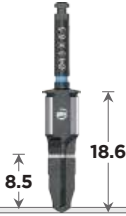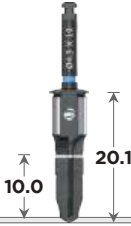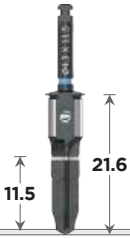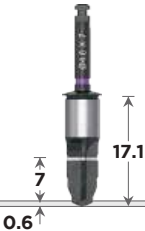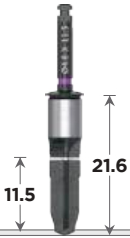Diameter		Length					
Scale 150% Unit : mm D Ø 4.3	Item No.						
		NSD43070	NSD43085	NSD43100	NSD43115	NSD43130	NSD43145
							
		NSD48070	NSD48085	NSD48100	NSD48115	NSD48130	

Socket Drill **Pylon**

Cuts the ridge of the extracted site by the fixture's diameter size, preventing drill slippage caused by remaining bone residue.
800 - 1200 rpm

Ø 3.5



ESD35

Ø 4.0



ESD40

Ø 4.5



ESD45

Ø 5.0



ESD50

Bone Profiler **Pylon**

Contours the bone around the coronal aspect to facilitate full seating of the abutment after removing the surgical guide.
800 - 1200 rpm

N



BPN4230

R



BPR5529

Extension Adapter

Extends the instrument length.

Scale 150%

Unit : mm

Handpiece (Latch)



SEAH40S



SEAH40L

***Sold Separately**

Ratchet



SEAR40S



SEAR40L

***Sold Separately**

Hex Driver 1.27

*Caution: Must be properly engaged and used at the recommended torque value.

Ratchet



ASDR127M



ASDR127L

Hand



ASD127M

Drill Extension

Extends the length of latch type tools.



DEH4314

Anchor Drill **Pylon**

Drills a hole on the bone for the anchor pin.



AD1413

Anchor Pin **Pylon**

Holds the surgical guide onto the bone.
Mostly used in an edentulous case.

***3 Anchor Pins Included**



AP15

Anchor Driver **Pylon**

Tightens down the anchor pin.



DMDR23

Ball Abutment Driver

Compatible with the ball abutment.

Ratchet



BADR24

Multi-Unit Abutment Driver

Compatible with Straight Multi-Unit abutment.

Ratchet



MUDR33

Fixture Driver **Pylon** A tool that delivers the fixture.

Platform			
Scale 150% Unit : mm N Narrow	Handpiece (Latch) GFDH20	Ratchet GFDR20S *Sold Separately GFDR20L	Ratchet (Stopper) GFDR20095 GFDR20110 GFDR20125 *Sold Separately
	Handpiece (Latch) GFDH25 *Caution: Not recommended on hard bone. Do NOT torque over 25N-cm and should not drill deeper than 6.5mm. Needs to be torqued with a ratchet driver as the final step.	Ratchet GFDR25S *Sold Separately GFDR25L Recommended Torque Value: 30 - 45N-cm	Ratchet (Stopper) GFDR25095 GFDR25110 GFDR25125 *Sold Separately Recommended Torque Value: 30 - 45N-cm
R Regular			

Torque Wrench

Tightens prosthetic components with the use of compatible ratchet-type drivers.
(Use at the recommended torque value of 20 - 40tN·cm)



TW40

Titanium Bowl




TB4762


Pylon Kit Guided Surgical Drilling Protocol

Implant Ø 5.0 x 10.0mm

Guide Sleeve Offset 9.5

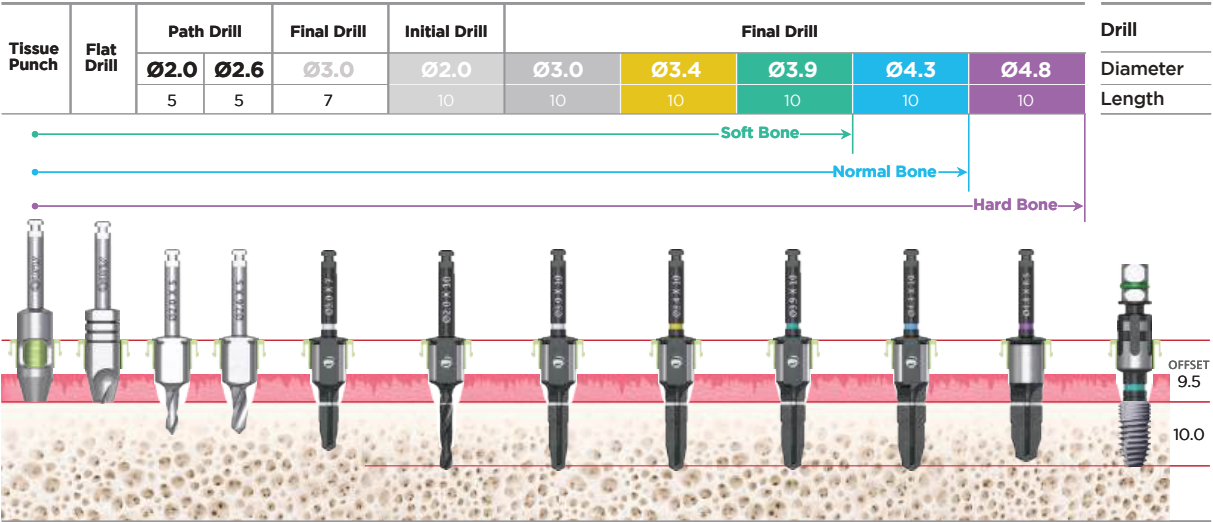


Guide Sleeve Offset 11.0

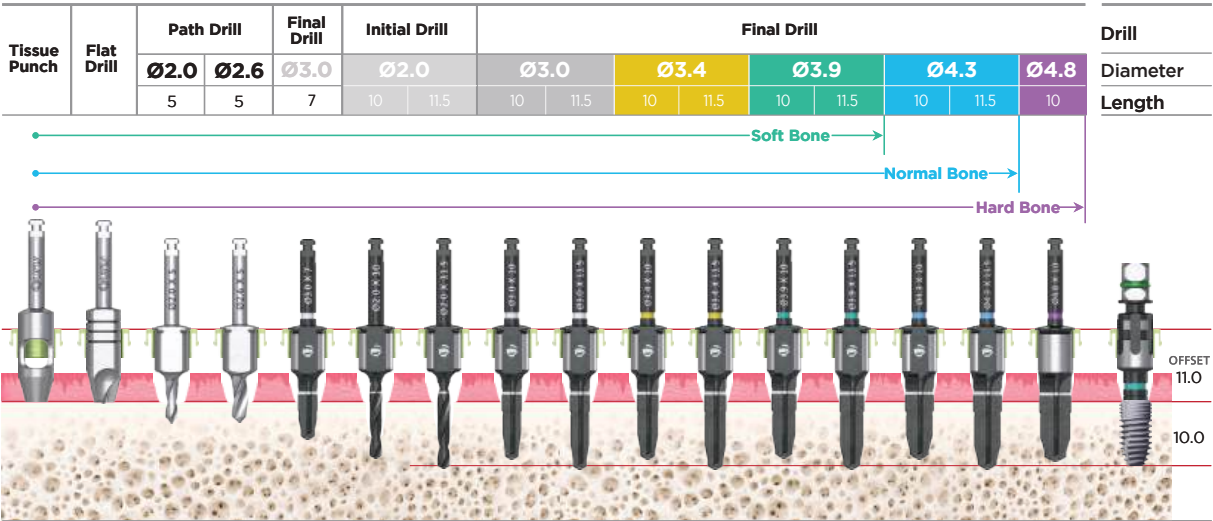

Guide Sleeve Offset 12.5

Implant Ø 5.0 x 10.0mm

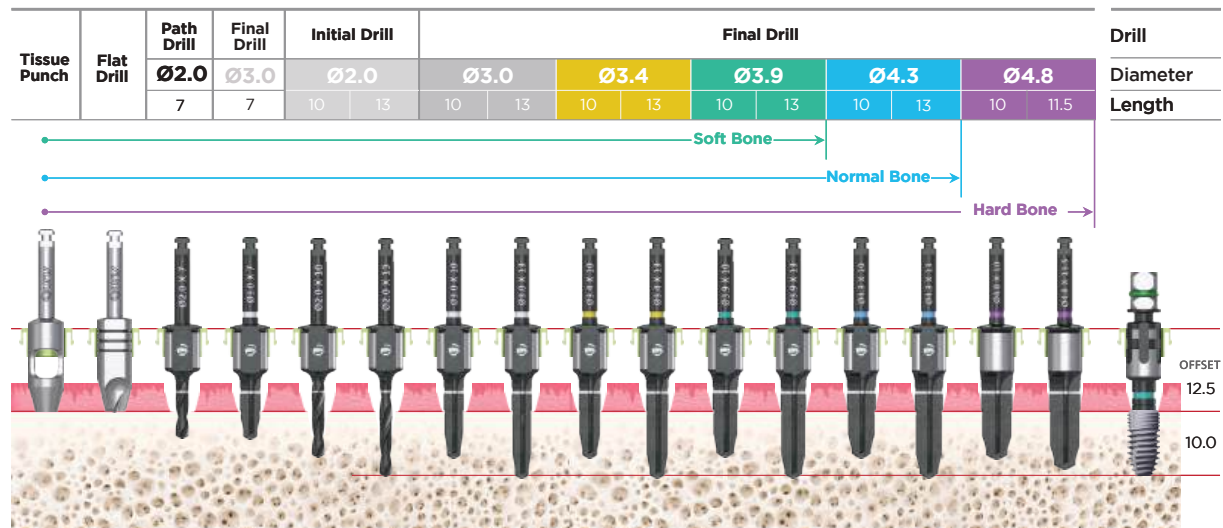

GUIDED PYLON PLUS KIT




PPK01



Pylon Plus Kit

Item No. PPK01

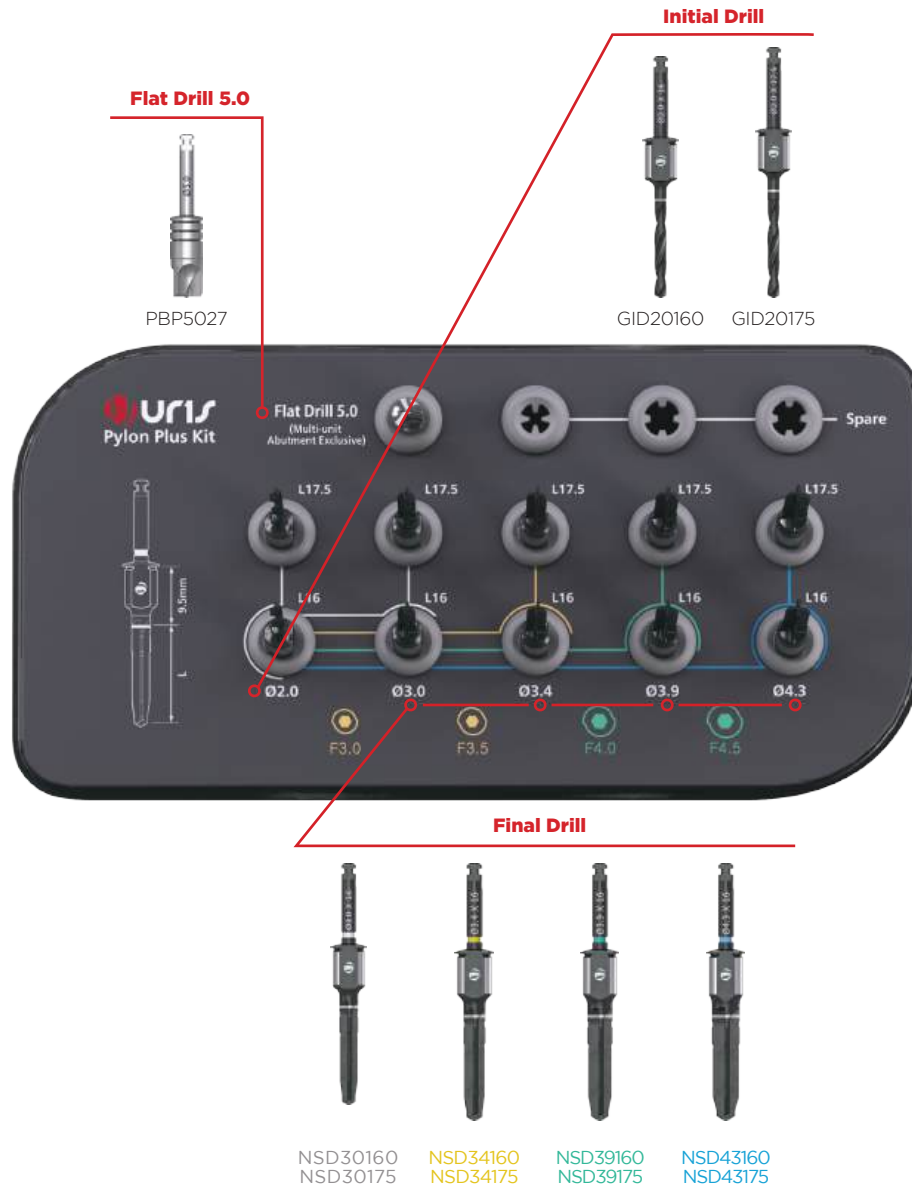


Flat Drill Pylon

Flattens out uneven bone and also removes remaining gingiva residue after tissue punching.

Scale 150%


Unit : mm

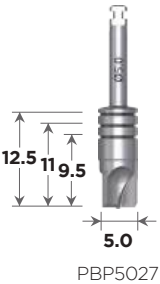

Initial Drill Pylon

Expands the drilling path by the chosen fixture's length.

Low speed: 50 - 100 rpm; within 5 seconds with high torque
High Speed: 800 - 1200 rpm (use irrigation)

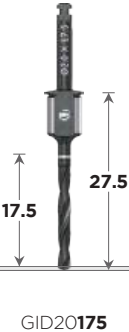Diameter	
 D Ø 2.0	 GID20160
	 GID20175
Item No.	

Final Drill Pylon

Used as the last drill to place a URIS fixture and is selected according to the fixture diameter and bone density.


Low speed: 50 - 100 rpm; within 5 seconds with high torque
High speed: 800 - 1200 rpm (use irrigation)

Scale 150%

Unit : mm

D Ø 3.0

Item No.

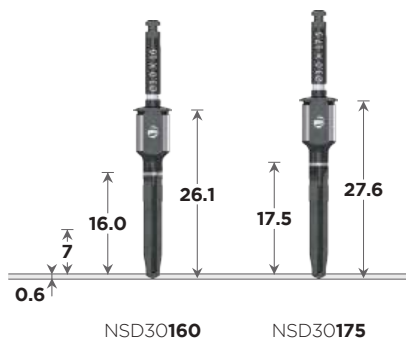


D Ø 4.3

Item No.



D Ø 3.4

Item No.

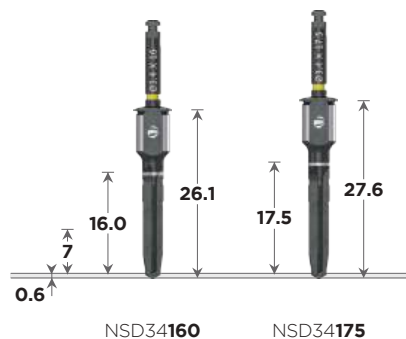

D Ø 3.9

Item No.



GUIDED


TAPERED PYLON KIT

TPK01



Tapered Pylon Kit

Item No. TPK01

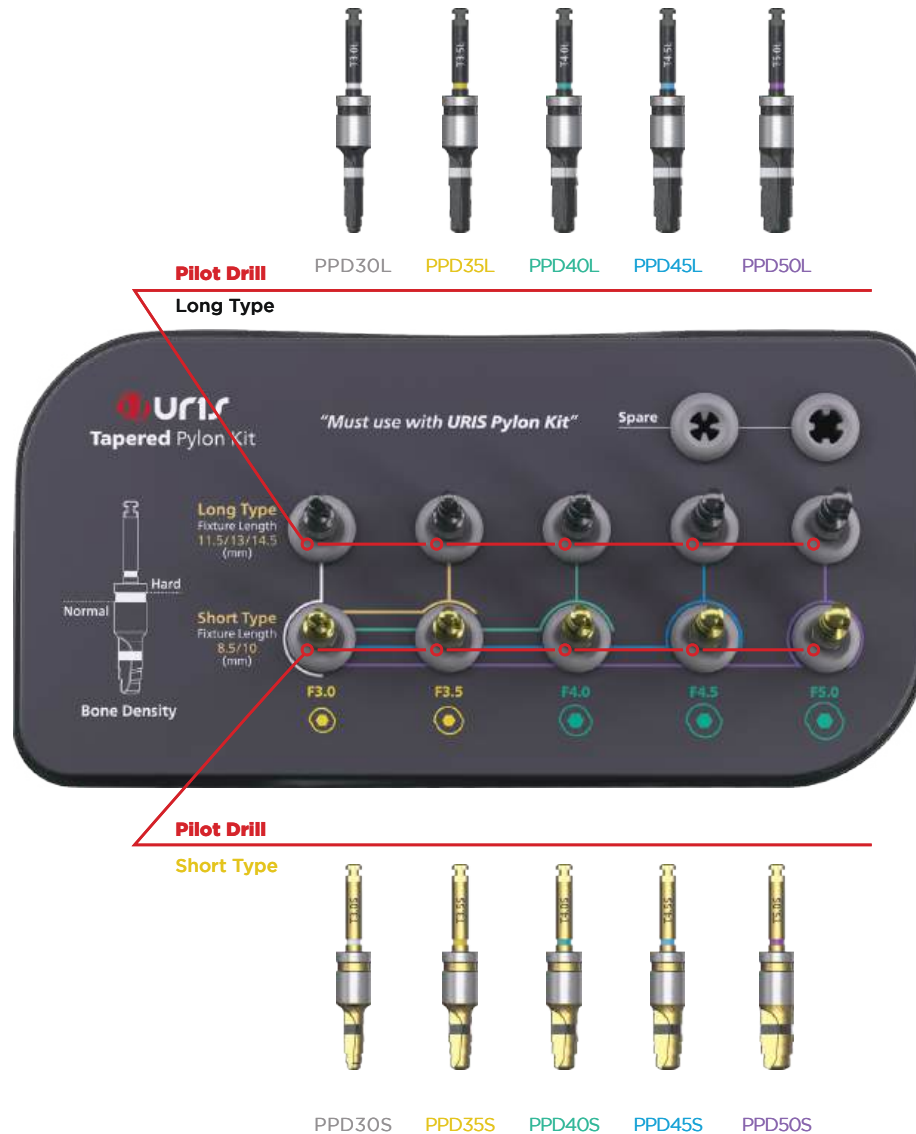

Pilot Drill Pylon

Used as the last drill to place a URIS tapered fixture and is selected according to the tapered fixture diameter and bone density.

Low Speed: 50 - 100 rpm; within 5 second with high torque
High Speed: 800 - 1200 rpm (use irrigation)

Length						
Scale 150% Unit : mm	Short Type	Ø 3.0	Ø 3.4	Ø 3.9	Ø 4.4	Ø 5.0
						
		Item No. PPD30S	PPD35S	PPD40S	PPD45S	PPD50S

Pilot Drill Pylon

Used as the last drill to place a URIS tapered fixture and is selected according to the tapered fixture diameter and bone density.

Low Speed: 50 - 100 rpm; within 5 second with high torque
High Speed: 800 - 1200 rpm (use irrigation)

Scale 150%
Unit : mm

Length					
Long Type	Ø 3.0	Ø 3.4	Ø 3.9	Ø 4.4	Ø 5.0
					
Item No.	PPD30L	PPD35L	PPD40L	PPD45L	PPD50L

GUIDED


NARROW PYLON KIT

UNP01



Narrow Pylon Kit

Item No. UNP01

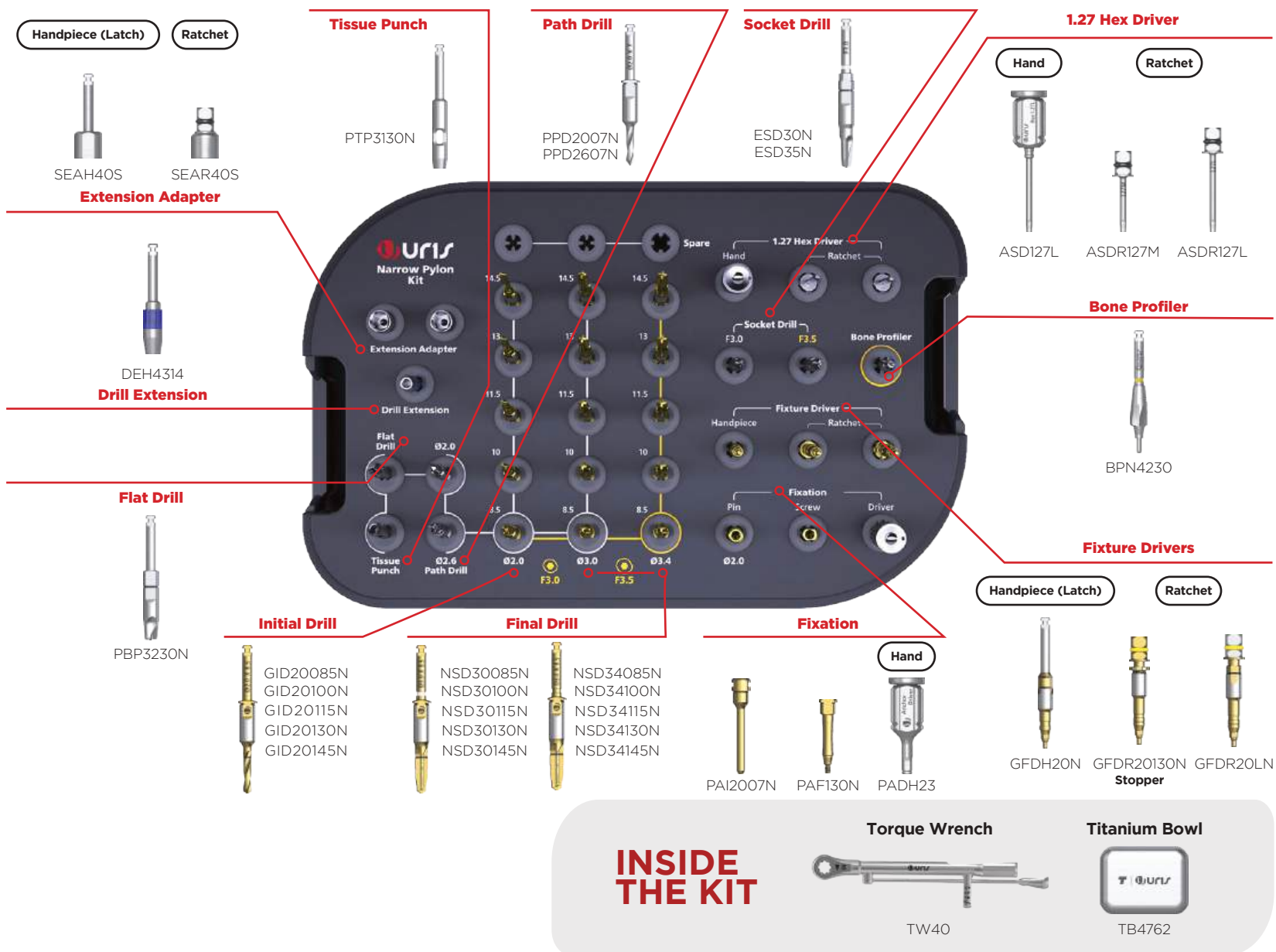

Tissue Punch Narrow Pylon

Removes soft tissue during a flap-less surgery.


Scale 150%

Unit : mm




Flat Drill Narrow Pylon

Flattens out uneven bone and also removes remaining gingiva residue after tissue punching.

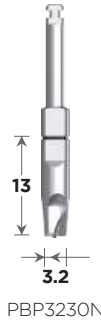

Path Drill Narrow Pylon

Increases the accuracy of the drilling path.

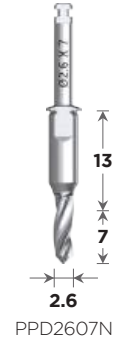

Tissue Punch/Path Drill/Initial Drill

Low speed: 50 - 100 rpm; within 5 seconds with high torque | High Speed: 800 - 1200 rpm (use irrigation)


Initial Drill Narrow Pylon

Expands the drilling path by the chosen fixture's length.

Low speed: 50 - 100 rpm; within 5 seconds with high torque
High Speed: 800 - 1200 rpm (use irrigation)

Implant Diameter

D Ø 2.0

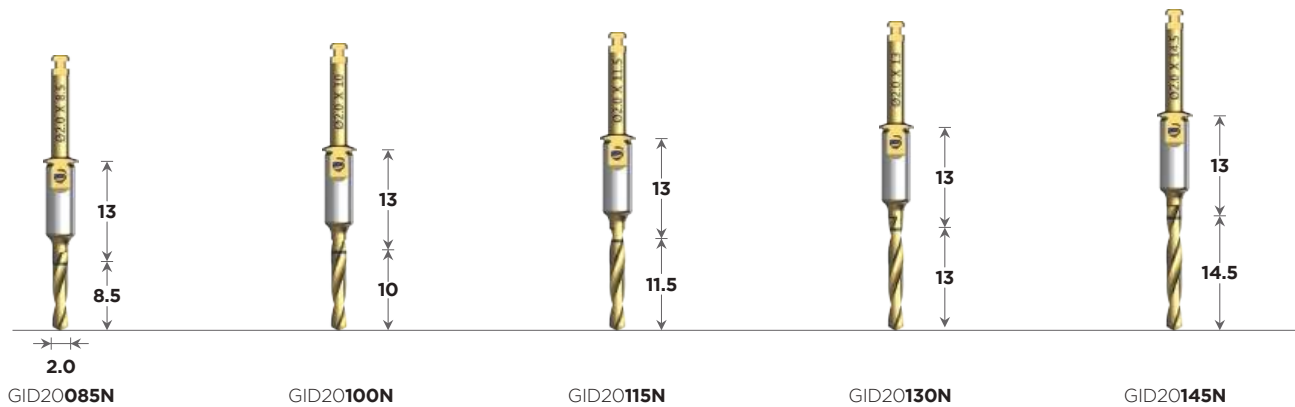

Item No.

Final Drill

Narrow Pylon

Used as the last drill to place a URIS fixture and is selected according to the fixture diameter and bone density.

Low speed: 50 - 100 rpm; within 5 seconds with high torque
High speed: 800 - 1200 rpm (use irrigation)

Implant Diameter

Length

Scale 150%

Unit : mm


D Ø 3.0

Item No.



D Ø 3.4

Item No.

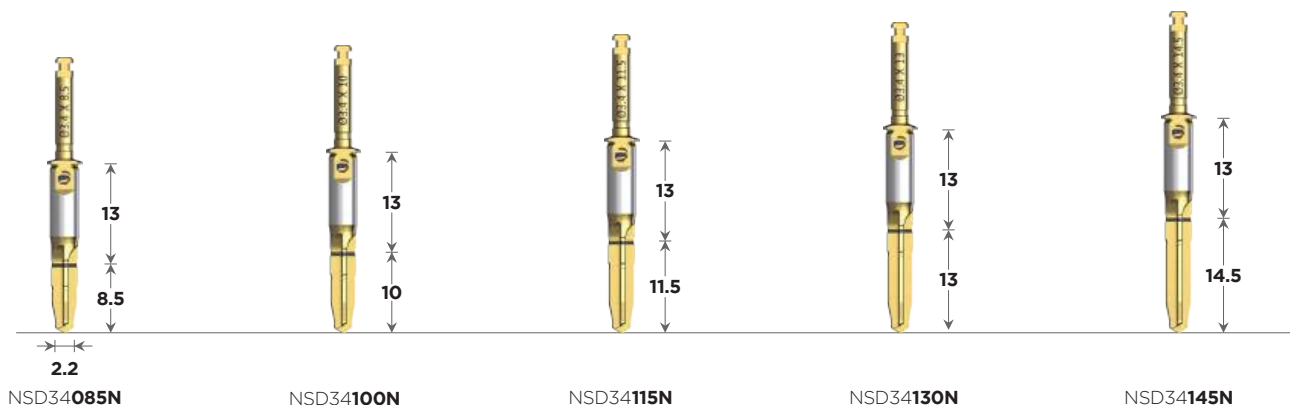

Socket Drill **Narrow Pylon**

Cuts the ridge of the extracted site by the fixture's diameter size, preventing drill slippage caused by remaining bone residue.
800 - 1200 rpm

Scale 150%

Unit : mm

Ø 3.0



ESD30N

Ø 3.5



ESD35N

Bone Profiler **Pylon**

Contours the bone around the coronal aspect to facilitate full seating of the abutment after removing the surgical guide.
800 - 1200 rpm



BPN4230

Fixture Driver **Narrow Pylon**

*Caution: Must torque after fully engaging into the fixture securely to avoid any breakages.
Handpiece: Up to 20N•cm
Ratchet recommended Torque Value: 30 - 45N•cm

Handpiece (Latch)

Ratchet (Stopper)

Ratchet



GFDH20N



GFDR20130N



GFDR20LN

Fixation Pin **Narrow Pylon**

Secures the surgical guide in place onto the osteotomy site after the 2.0 drill on edentulous cases to deliver the fixtures.



PAI2007N

Fixation Screw **Narrow Pylon**

Holds the surgical guide after the fixture delivery.

Scale 150%

Unit : mm



PAF130N

Anchor Driver **Pylon**

Tightens down the anchor pin.

Hand



PADH23

Extension Adapter

Extends the instrument length.

Handpiece (Latch)



SEAH40S

Ratchet



SEAR40S

Drill Extension

Extends the length of latch type tools.



DEH4314

Hex Driver 1.27

*Caution: Must be properly engaged and used at the recommended torque value.

Scale 150%

Unit : mm

Ratchet

Hand



ASDR127M



ASDR127L



ASD127L

Torque Wrench

Tightens prosthetic components with the use of compatible ratchet-type drivers.
(Use at the recommended torque value of 20 - 40N·cm)



TW40



Titanium Bowl



TB4762




GUIDED

PYLON ANCHOR KIT

PAK01


Pylon Anchor Kit Item No. PAK01

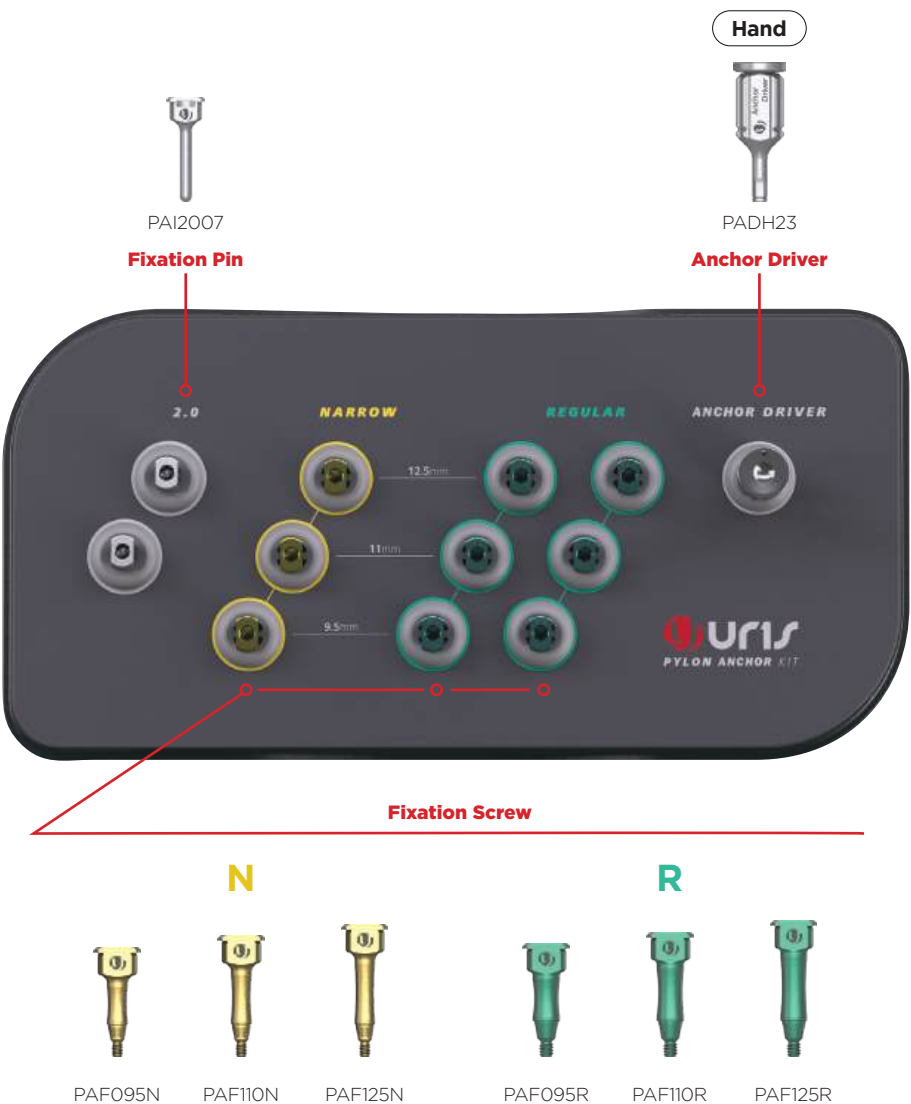

Fixation Pin Pylon

Secures the surgical guide in place onto the osteotomy site after the 2.0 drill on edentulous cases to deliver the fixtures.

***2 Fixation Pins Included**

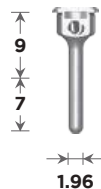

PAI2007

Anchor Driver Pylon

Tightens down the anchor pin.

Hand




PADH23




Fixation Screw Pylon

Holds the surgical guide after the fixture delivery.

N

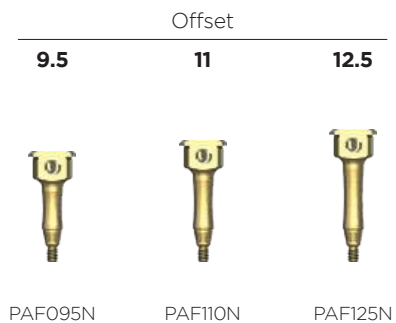

R

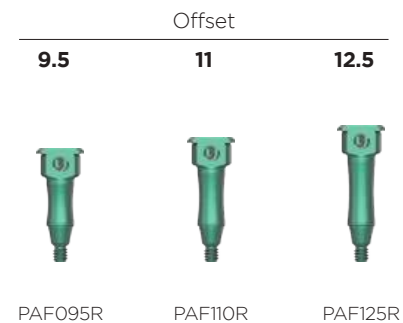

GUIDED


PYLON CRESTAL SINUS KIT

PSK01



Pylon Crestal Sinus Kit

Item No. PSK01

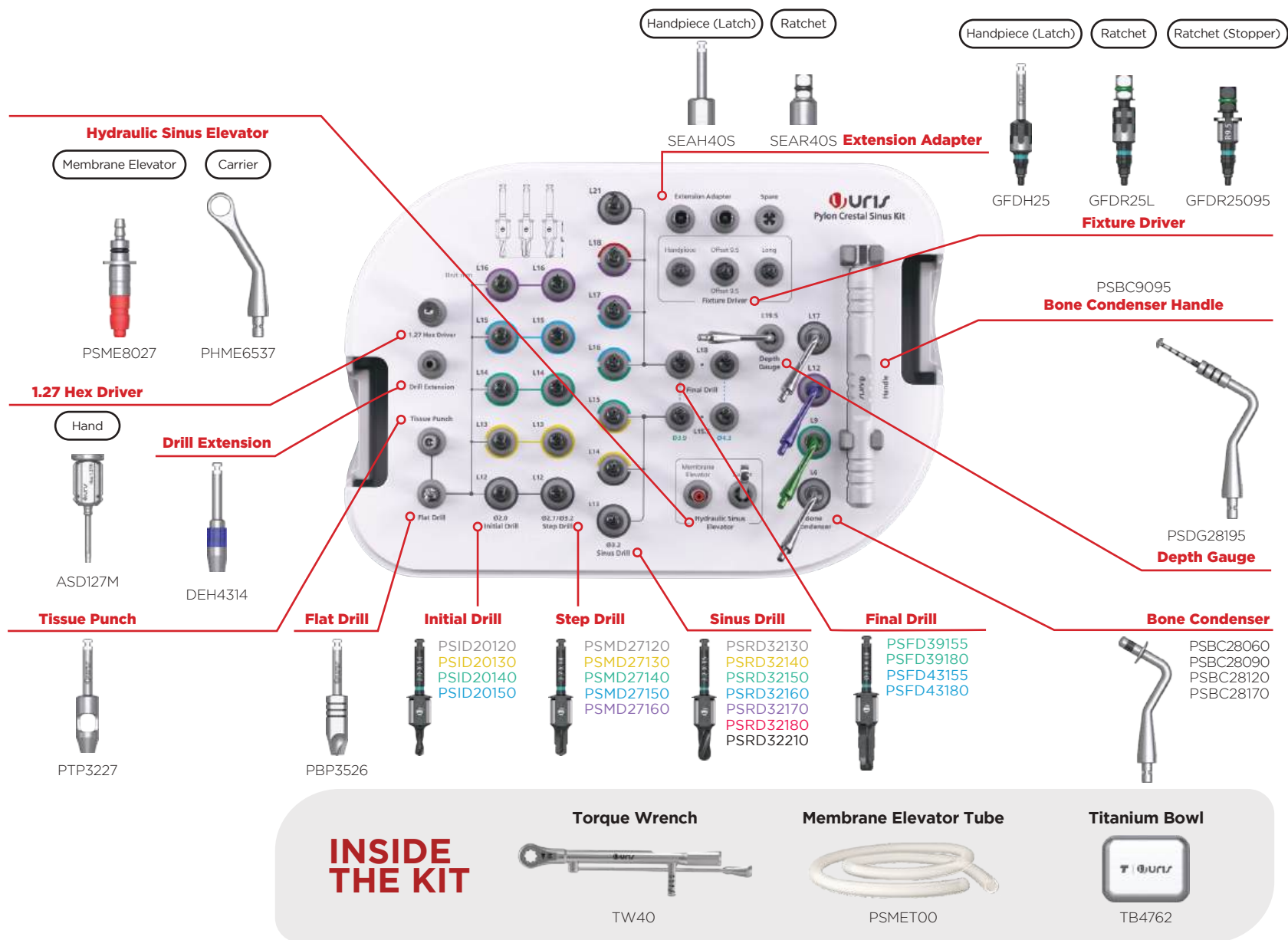

Initial Drill Pylon Sinus

Enlarges the osteotomy site. Do NOT drill deeper than 0.5 - 1.0mm directly below the maxillary sinus floor.

Low speed: 50 - 100 rpm; within 5 seconds with high torque
High speed: 800 - 1200 rpm (use irrigation)

Diameter					
Scale 150% Unit : mm					
		D Ø 2.0			
Item No.					
		PSID20120	PSID20130	PSID20140	PSID20150
					PSID20160

Step Drill Pylon Sinus

Used after D2.0 initial drill to enlarge the osteotomy site. DO NOT drill deeper than 0.5 - 1.0mm directly below the maxillary sinus floor.

Low speed: 50 - 100 rpm; within 5 seconds with high torque
High speed: 800 - 1200 rpm (use irrigation)

Diameter					
Scale 150% Unit : mm					
		D Ø 3.0			
Item No.					
		PSMD20120	PSMD20130	PSMD20140	PSMD20150
					PSMD20160

Sinus Drill

Pylon Sinus

Used to approach the sinus membrane and should drill 1 - 2mm deeper than the step drill.

Low speed: 50 - 100 rpm; within 5 seconds
with high torque (*Drill while applying pressure)
High speed: PROHIBITED

Scale 150%

Unit : mm

D Ø 3.2

Item No.

Diameter



PSRD32130



PSRD32140



PSRD32150



PSRD32160



PSRD32170



PSRD32180



PSRD32210

Final Drill

Pylon Sinus

Used with a surgical guide after bone grafting.

Low speed: 50 - 100 rpm; within 5 seconds with high torque

D Ø 3.9

Item No.

Diameter



PSFD39155

PSFD39180

D Ø 4.3

Diameter



PSFD43155

PSFD43180

Hex Driver 1.27

*Caution: Must be properly engaged and used at the recommended torque value.

Scale 150%

Unit : mm

Hand



ASD127M

Drill Extension

Extends the length of latch type tools.



DEH4314

Tissue Punch


Removes soft tissue during a flap-less surgery.
800 - 1200 rpm



PTP3227

Flat Drill Pylon

Flattens out uneven bone and also removes remaining gingiva residue after tissue punching.

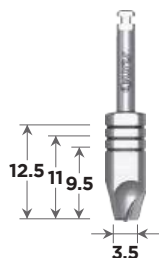

PBP3526

Fixture Driver

Pylon

A tool that delivers the fixture.

Scale 150%

Unit : mm

R
Regular

Platform

Handpiece (Latch)



GFDH25

*Caution: Not recommended on hard bone.
Do NOT torque over 25N-cm and should not
drill deeper than 6.5mm. Needs to be torqued
with a ratchet driver as the final step.

Ratchet



GFDR25L

Recommended Torque Value: 30 - 45N-cm

Ratchet (Stopper)



GFDR25095

Recommended Torque Value: 30 - 45N-cm

Membrane Elevator Carrier

Pylon Sinus

Delivers the membrane elevator.
Needs to be engaged onto a handle.

Scale 150%

Unit : mm



PHME6537

Membrane Elevator

Pylon Sinus

Lifts sinus membrane with hydraulic pressure.
Needs to be engaged onto a carrier.



PSME8027

Extension Adapter

Extends the instrument length.

Handpiece (Latch)

Ratchet



SEAH40S



SEAR40S

Membrane Elevator Tube

A transparent silicon tube.
Outer D4.0/Inner D2.0/Length 300mm | Autoclave before use | Single use only



PSMET00

Titanium Bowl




TB4762

Bone Condenser Pylon Sinus

Insert bone graft materials into interior wall once the wall is completely opened. Used after removing the surgical guide.

Scale 150%

Unit : mm

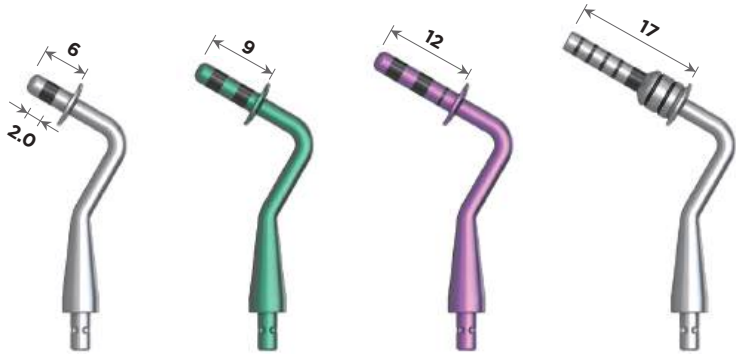

PSBC28060

PSBC28090

PSBC28120

PSBC28170

Bone Condenser Handle Pylon Sinus



PSBC9095

Depth Gauge Pylon Sinus

Measures depth with measurement markings at the tip.



PSDG28195

Torque Wrench

Tightens prosthetic components with the use of compatible ratchet-type drivers.
(Use at the recommended torque value of 20 - 40N.cm)




TW40


Pylon Crestal Sinus Kit Guided Surgical Drilling Protocol

Implant Ø5.0



Guide Sleeve Offset 9.5

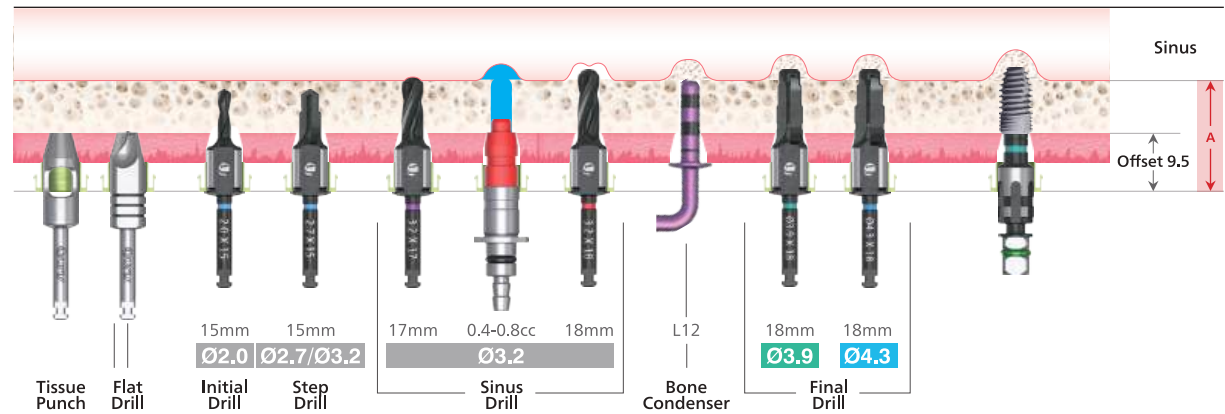

Guide Sleeve Offset 11.0

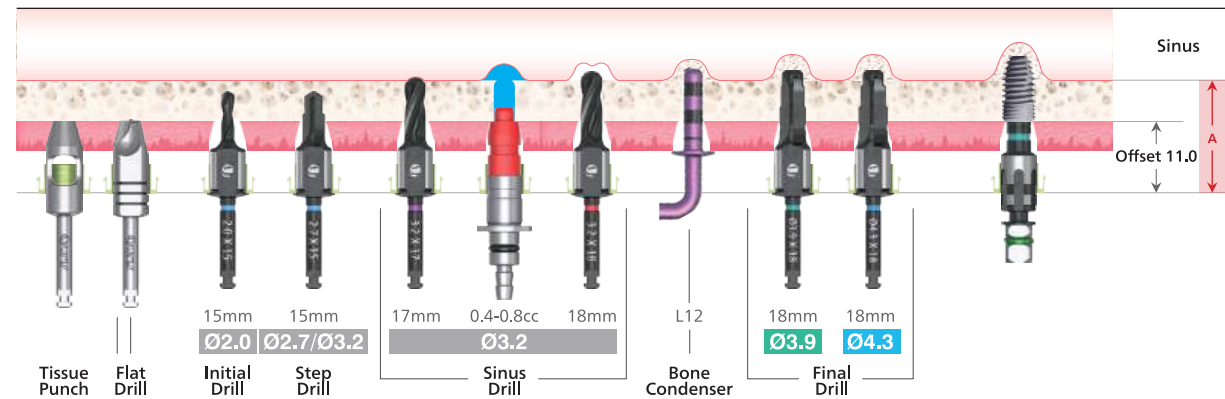

Guide Sleeve Offset 12.5

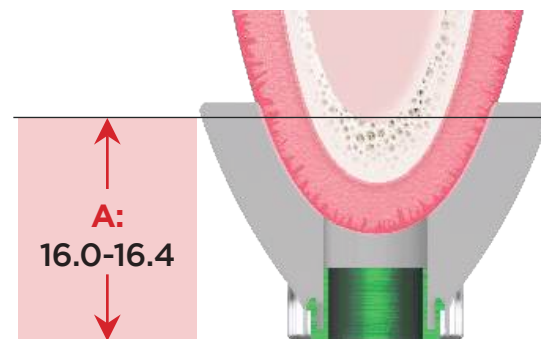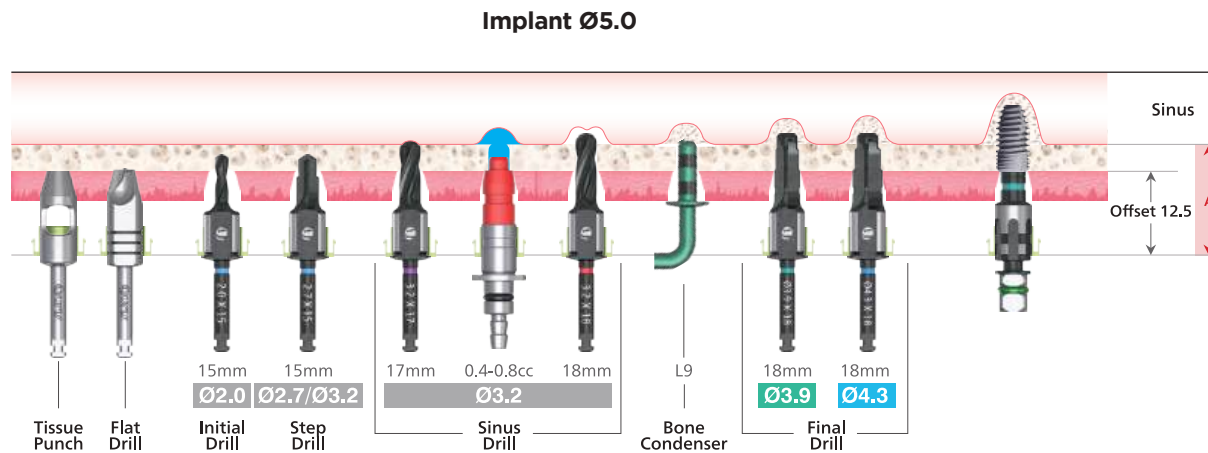


PROSTHETIC DRIVER KIT

PDK01



Prosthetic Driver Kit Item No. PDK01

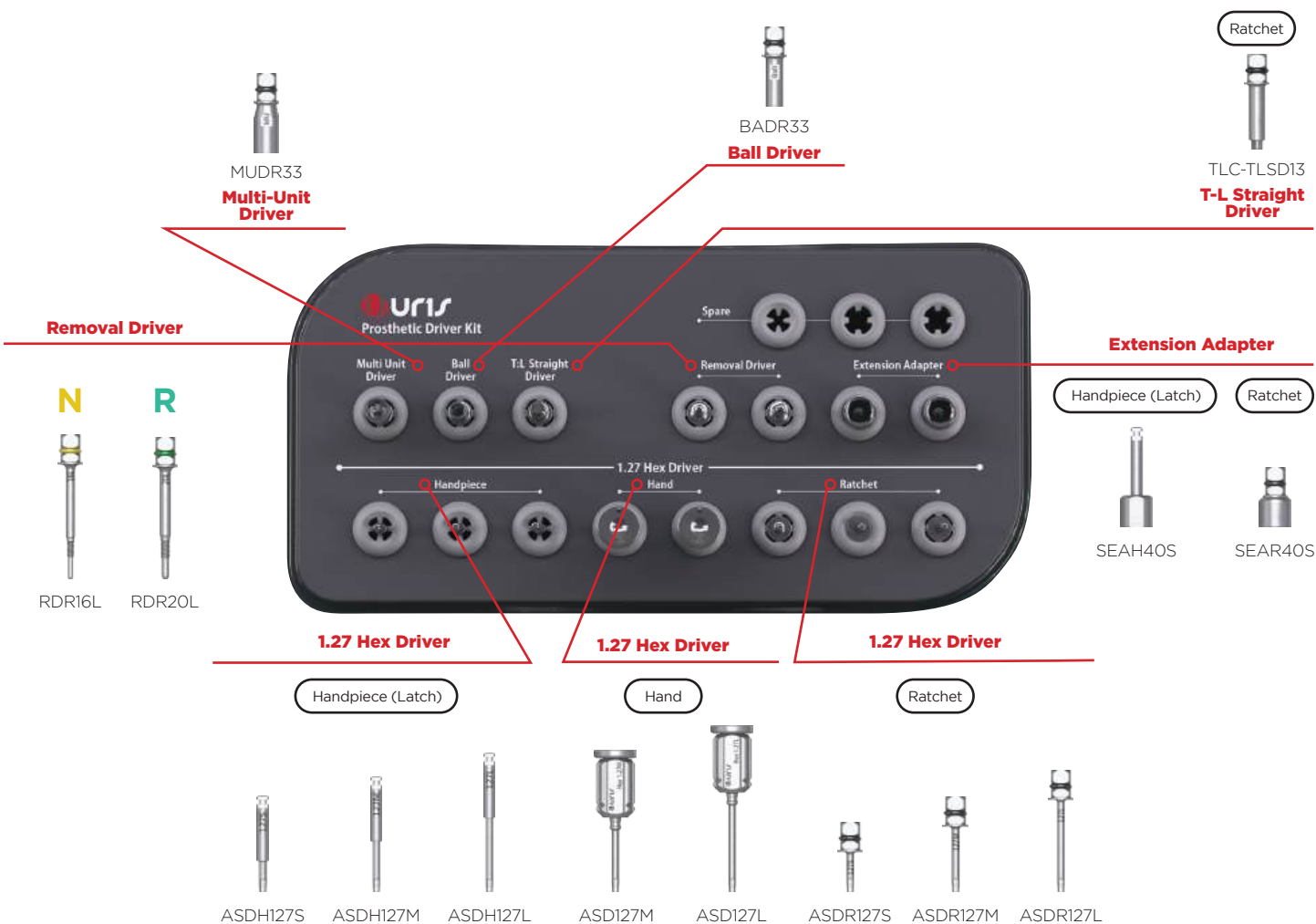

**INSIDE
THE KIT**

Torque Wrench



TW40

Hex Driver 1.27

*Caution: Must be properly engaged and used at the recommended torque value.

Scale 150%

Unit : mm

Hand



ASD127M



ASD127L

Handpiece (Latch)



ASDH127S



ASDH127M



ASDH127L



ASDH127LL
*Sold Separately

Ratchet



ASDR127S



ASDR127M



ASDR127L



ASDR127LL
*Sold Separately

Ball Abutment Driver

Compatible with Ball abutment.

Ratchet



BADR24

Multi-Unit Abutment Driver

Compatible with Straight Multi-Unit abutment.

Ratchet



MUDR33

T-L Straight Abutment Driver

Compatible with T-L abutments.

Ratchet



TLC-TLSD13

Scale 150%

Unit : mm

Torque Wrench

Tightens prosthetic components with the use of compatible ratchet-type drivers.
(Use at the recommended torque value of 20 - 40N·cm)



TW40

Removal Driver

Used to remove the abutment off the fixture.
Needs to be inserted upright and tightened clockwise.

N

R



RDR16S

*Sold Separately



RDR16L



RDR20S

*Sold Separately



RDR20L

Extension Adapter

Extends the instrument length.

Handpiece (Latch)

Ratchet



SEAH40S




SEAR40S

ASC KIT

ASC01




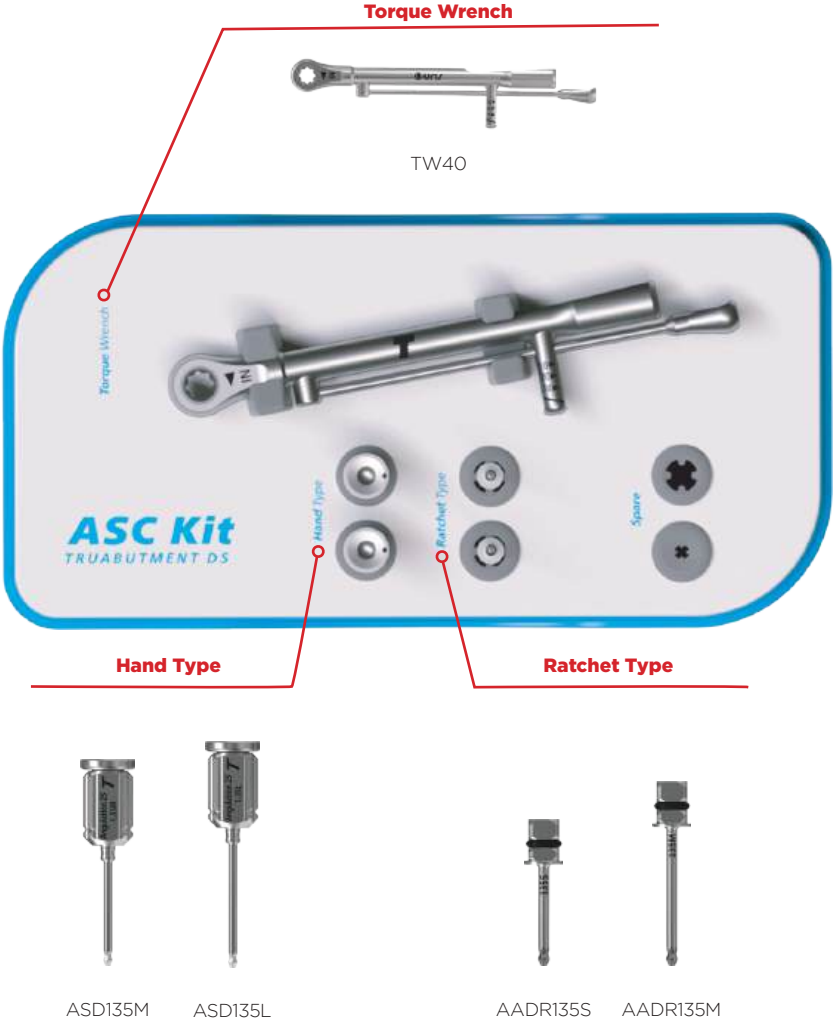

(ASC) Angulated Screw Channel Driver

*Caution: Must be properly engaged and torqued at the recommended torque value of 20 - 25N·cm within 25°.

Hand

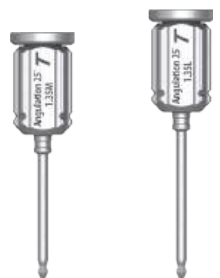

ASD135M

ASD135L

Ratchet



AADR 135S

AADR 135M



AADR135L

***Sold Separately**

Handpiece (Latch)



AADH135S



AADH135M



AADH135L

***Sold Separately**


Torque Wrench

Tightens prosthetic components with the use of compatible ratchet-type drivers.
(Use at the recommended torque value of 20 - 40N·cm)



TW40

Fixture Packaging

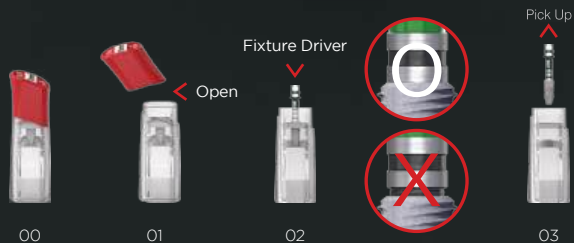


Cover Screw



Cover Cap
Cover Screw
Upper Cap

Abutment Packaging

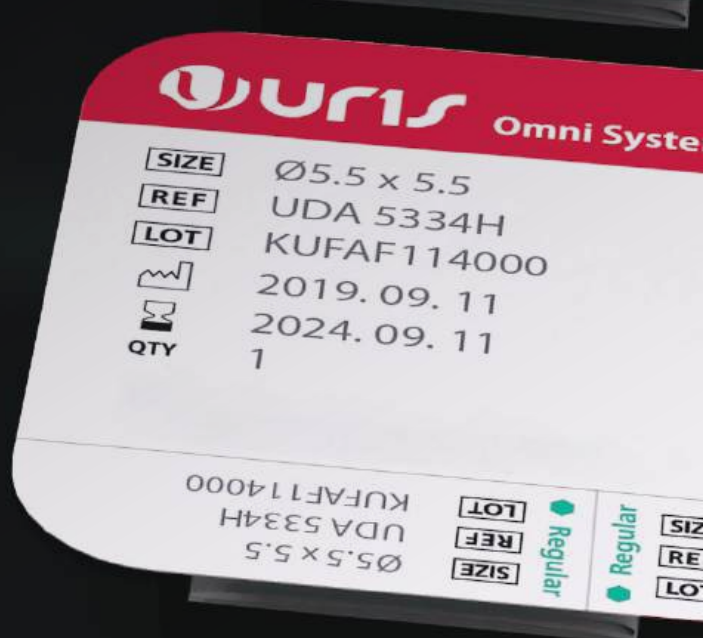



Instructions for use

SCAN ME FOR IFU



TruAbutment Korea Co., Ltd.

#214, 215, 216, 104-dong, Seokcheon-ro 397,

Ojeoung-gu, Bucheon-si, Gyeonggi-do, Korea

Phone: +82 (32) 678 4688

Fax: +82 (32) 624 4689

www.truabutment.com

www.urisimplants.com



Ultimate Restorative Implant Solution

URISIMPLANTS.COM